

R.O.C

Banhu Industrial District, Huang Jiang, Dong Guan City, Guan Dong Province, P.R.C (PLANT)

86-769-83361381

THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	Protection ¹	Respirator	Cautions and Limitations ²
		MF-01	
84A-7697	Surgical/N95	X	ABCJMNO

1. Protection

N95 - Particulate Filter (95% filter efficiency level)
 Effective against particulate aerosols free of oil;
 time use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to user instructions, and/or maintenance manuals for information on use and maintenance of these respirators.

S - Special or critical User Instructions and/or specific use limitations apply. Refer to User Instructions before donning.


This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings, as a Surgical N95 Respirator

conforming to recognized standards for biocompatibility, flammability, and fluid resistance.

Note :

- 1 - Not for use in atmospheres containing less than 19.5% oxygen.
- 2 - Not for use in atmospheres immediately dangerous to life or health.
- 3 - Do not exceed maximum use concentrations established by regulatory standards.
- 4- Failure to properly use and maintain this product could result in injury or death.
- 5 - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- 6 - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- 7 - Refer to user instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- 8 - Special or critical User Instructions and/or specific use limitations apply. Refer to User Instructions before donning.**
- 9- Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
- 10-Storage conditions: place in dried and ventilated place, do not make the sun bask directly.
- 11-Please pay attention to the validity of the bags (boxes) before use, over the period of validity is strictly prohibited.
- 12-The respirators are disposable. Can not reuse.
- 13-The product can't contact with the flame and the extremely hot object in order to avoid burning.
- 14-Please check the mask whether surface is damaged, ties fracture and attach insecurely and nose clip (nose piece) exposes outside before opening the package. If have, do not use.
- 15-Please refer to the instruction for use when wear.
- 16- If you cannot achieve a proper fit, DO NOT enter the contaminated area. See your supervisor or call 0086-0769-83361381 Dong Guan Master & Frank Enterprise Co., Ltd.

Donning & Fitting Instruction

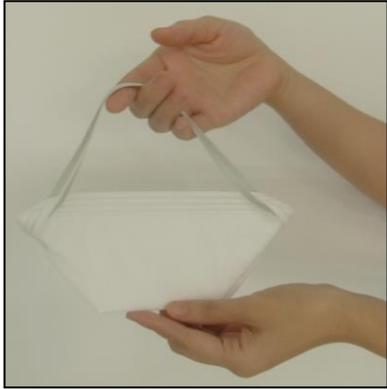

1. Adjust the elastic straps for suitable size before placing the respirator on the face.



2. Cup the respirator in your hand, with nose piece towards your fingers.



3. Position the respirator at mouth and nose. Pull the elastic strap near your hand over your head resting it at the back of your head. Pull the strap near your finger over your head and position it at the top of your head and adjust the elastic straps.



4. Shape the nose piece suitably with your forefingers and middle fingers in order to keep the respirator to cling your face.



5. Use two hands perform a user seal check prior to each wearing. To check fit, place both hands completely over the respirator and exhale. If air leaks around nose, readjust the nosepiece as described in step 4. If air leaks at the respirator edges, work the straps back along the sides of your head.