


The NIOSH approval label will be printed on the box.

**N95-5150 NIOSH  
N95 Particulate Respirator**

**User's Instruction Manual**

**Dromex CC  
Unit 1, Blasé Road New Germany  
South Africa, 3610  
+27 (31) 713 1960**

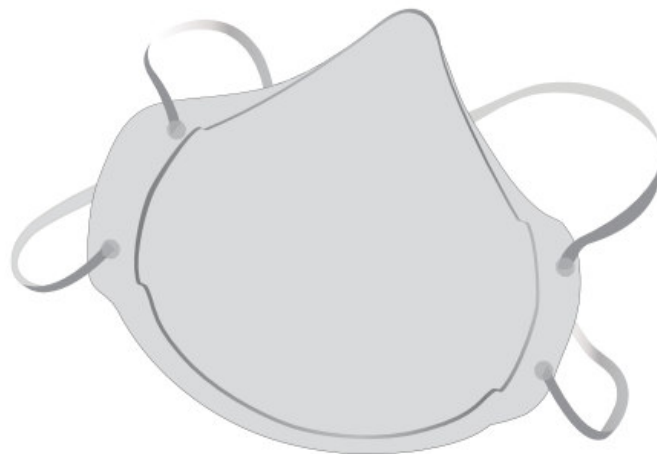

**N95: 95% filtration efficiency against solid and liquid aerosol that do not contain oil.**

**MEETS CDC GUIDLINES FOR M.TUBERCULOSIS EXPOSURE CONTROL**

**ATTENTION: READ APPLICATIONS, LIMITATIONS AND FITTING INSTRUCTIONS ON BOX BEFORE RESPIRATOR IS WORN.**

**CONTENTS: 20 RESPIRATORS**

# **APPLICATIONS**

## **USE AGAINST:**

NIOSH approved for protection against airborne particulate aerosols free of oils.

These masks are suitable for use in occupational areas such as mining, agriculture, construction, general material handling and medical industry applications. They are also suitable as one part of an infection-control strategy.

## **DO NOT USE AGAINST:**

1. Vapors, paint spray, oil-based aerosols, asbestos and sand blasting.
2. Any contaminant whose concentration is unknown, or immediately dangerous to life.
3. Concentrations of Dust and Mists which exceed the maximum use concentration or 10 times contaminant's permissible exposure limit established by OSHA.

## **WARNING:**

1. These respirators do not supply oxygen. Use only in well ventilated areas containing sufficient oxygen to support life.
2. Leave contaminated area immediately if the respirator becomes damaged and proper fit cannot be obtained or breathing becomes difficult. Dispose and change the respirator.

## **IMPORTANT LIMITATIONS:**

1. All instructions and warnings have to be followed. Failure to do so will reduce the efficiency of respirator and may cause sickness or death.
2. Training on the proper use of respirator must be provided by the employer in accordance with the safety and health standard.
3. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
4. Do not alter or modify the respirator.
5. Not suitable for persons with beard or other facial hair which may reduce the efficiency of respirator.
6. The use and reuse of this respirator is subject to consideration of hygiene damage and increased breathing resistance. For dirty workplaces, service time should not be extended beyond 8 total hours of use.

7. Store the respirators in cool place. Avoid direct exposure of sunlight during storage as the filter media in this respirator can be degraded by extended exposure to direct sunlight.
8. S-Special or critical User Instructions: The model N95-5150 filter facepiece respirator has been manufactured by JINFUYU INDUSTRIAL CO., LTD. for DROMEX under TC-84A-7719.

## Instructions for the N95-5150 Maintenance free respirator

### Fitting Instructions

The respirator user must be fit tested according to OSHA regulation 29 CFR 1910.134.

1. Hold the respirator in one hand with the nose area at your fingertips, allowing the headbands to hang below your hand.
2. Press the respirator against your face with the nose foam on the bridge of your nose.
3. Place the bottom band over your head and position it below your ears. Pull the top strap high on the back of your head.
4. Test the fit. Cup both hands over the respirator and exhale vigorously. If any leakage is detected, readjust the position of the respirator according to Steps 2, 3 & 4 to achieve a better fit.

If you CANNOT achieve a proper fit, DO NOT enter the contaminated area.  
See your supervisor for proper fitting respiratory protection.

