

USER'S INSTRUCTIONS

MC Medical, Inc.
N95 Respirator

MW09004NM (Medium)
MW09004NS (Small)
MW9004NMO (Medium)
MW9004NSO (Small)
MW9004NMB (Medium)
MW9004NSB (Small)
MW9004NMG (Medium)
MW9004NSG (Small)

MC Medical, Inc.
Shinagawa East One Tower Floor 12th, 2-16-1, Konan,
Minato-ku, Tokyo, 108-0075 Japan
TEL : +81-3-5715-2818

TN0206PL3U1a
[Rev. a] 02/08/2021

USER'S INSTRUCTIONS

WARNING!

Although this respirator is useful from particulate aerosols free of oil, the risk for sickness or infection is not completely eliminated. Failure to read and follow all warnings and instruction will reduce respirator performance and may result in sickness and / or death.

[USE INSTRUCTION]

This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.

Use for protection against particulate aerosols free of oil.

Do not use for gasses, vapors asbestos, paint spray or any oil-based aerosols, or in environments.

Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

Leave the use environment immediately in the case of dizziness, nausea, breathing difficulty or other symptoms of physical distress.

Don a new respirator if a proper fit cannot be obtained, if breathing becomes difficult, or the respirator becomes damaged.

Only for single use.


Store unused respirators in a clean, cool, dark, dry place.

S-[SPECIAL OR CRITICAL USER'S INSTRUCTIONS]


These MW09004NM and MW09004NS, MW9004NMO, MW9004NSO, MW9004NMB, MW9004NSB, MW9004NMG, MW9004NSG N95 filtering face piece respirators" have been manufactured by Approval Holder San-M Package Co., Ltd. for MC Medical, Inc. under TC-84A-8542.

Fitting Instruction


* Donning procedure

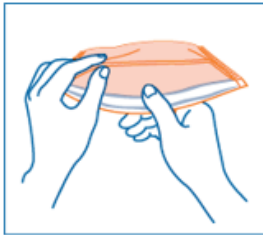

1. Please inspect the respirator prior to each use to assure there is no damage, and then separate mask to open.

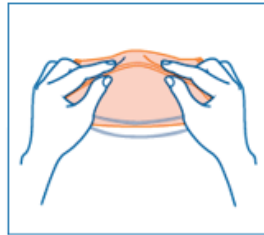

2. Bend the nosepiece to conform respirator to your nose.

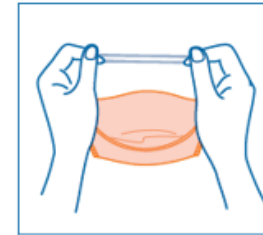

3. Hold the center of headbands separately.



4. Hang the lower part on your chin and bring both bands over your head.



5. Place the upper band at your crown and the lower band at your neck.



6. Form the nosepiece tightly to your face by using the fingers of both hands, form the nose piece to the nose.

※ To check respirator fit: Cup both hands around the respirator, being careful not to disturb the position. Breathe out, if no breath escapes from around the nose, cheeks or edges of the respirator the respirator is fitted properly. If air leaks around the facepiece, reposition the straps and try again. Follow these fitting instructions during each use to achieve a proper fit.

* Doffing procedure



1. Grab both elastic band with both hands, and pull it backwards to remove it from the head.