

FLAT FOLD DISPOSABLE RESPIRATOR

PROTECTION

N95-Particulate Filter (95% filter efficiency level). Effective against particulate aerosols free of oil. Time use restrictions may apply.



FITTING INSTRUCTIONS

Follow these instructions each time respirator is worn.

1. Open the respirator by pushing slightly on both sides.
2. Place the mask under the chin, the nose-bridge upward and outside.
3. Place the lower headband around the neck, below the ears. Holding the mask against the face with one hand, place the top headband above the ears, around the crown of the head.
4. Mold the nose area to the shape of the nose by pinching the nose-bridge from the top to the bottom.

FIT-CHECK THE FACE-SEAL AS FOLLOWS:

1. Place both hands over the respirator - do not disturb the position of the respirator;
2. Inhale sharply (a negative pressure should be felt inside the respirator);
3. If you detect air-leaks, readjust the head straps and/or the nose-piece.

WARNING: If you cannot achieve satisfactory fit, do NOT enter the contaminated area. Consult your Super visor.

USER INSTRUCTIONS:

1. Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the whole duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury or death.
2. Before Occupational Use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirement. **In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing, and applicable OSHA substance specific standards.**
3. Discard the respirator and replace with a new one if:
 - Excessive clogging of the respirator causes breathing difficulty.
 - The respirator becomes damaged.

4. Visually inspect the filtering facepiece prior to each use to assure that there are no holes in the breathing zone, especially around the weld attachment of the headbands to insure no damage has occurred. Discard respirator if any visual signs of holes or damage. **Please inspect the headbands to insure that they are firmly attached to the respirator.**
5. Leave the contaminated area if dizziness, irritation or other distress occurs.

USE LIMITATIONS:

1. Do NOT use the respirator or enter or stay in a contaminated area under the following circumstances:
 - Atmosphere contains oil aerosols
 - For protection against gases or vapors
 - Contaminants or their concentrations are unknown or immediately dangerous to life or health
 - Concentrations or contaminants exceed maximum use concentrations in applicable OSHA standards or applicable government regulations or 10 times the PEL (Permissible Exposure Limit), whichever is lower
 - For sandblasting, asbestos
2. Do NOT modify or misuse the respirator.
3. Do NOT use the respirator with beards or other facial hair that interferes with direct contact between the face and the edge of the respirator, or any other conditions that may prevent a good face-seal.


WARNING: This respirator helps protect against certain contaminants. It may not eliminate the risk of contracting disease or infection. Misuse may result in sickness or death.

Materials which may come into contact with the wearer's skin could cause allergic reactions to susceptible individual. FOR MORE INFORMATION or assistance on Honeywell products, contact your local Honey well distributor.

STORAGE CONDITIONS:

Temperature: Between -20°C to 38°C.

Humidity: Less than 80%

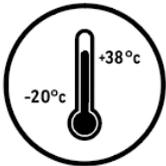

Contact Information for this product: Indrajit.Ghorpade@honeywell.com