

N95 Filtering Facepiece

Foldable

Single-Use

Model: M8004

Very close facial fit

>95% filtration efficiency against solid and liquid aerosols free of oil

### Description

The Particulate Respirator M8004 is designed to help provide respiratory protection for the user. This product has a filter efficiency of 95% or greater against solid and liquid aerosols free of oil. This product contains no components made from natural rubber latex.

### Use Instructions


1. Suitable for protecting the mouth, nose, and chin to airborne particles and block liquid aerosols.
2. For adults use only.
3. Do not wash.
4. Store in a cool, dry, clean place away from fire and contamination.
5. Expiration date: 2 years after production.

### 4-Step check before use

1. Do not use if packaging has been opened or damaged.
2. Check if elastic bands are in good condition and are not damaged.
3. Check that the nose clip is not broken.
4. Check there are no holes or damage in the breathing zone.

### Wearing instructions

1. Unfold the respirator and hold it with both hands, with the nose clip facing up.
2. Place elastic bands around the neck and head respectively. Put the respirator against your face covering both nose and mouth.
3. Adjust the nose clip cover bridge of the nose using two fingers to press down until achieving a close fit.
4. Perform a fit check according to instructions.

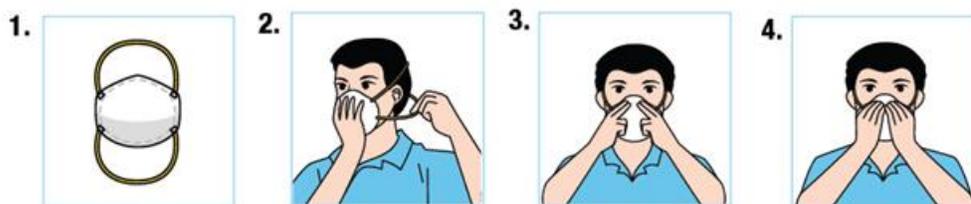

### Fit check

1. Place both hands over respirator.
2. Take a deep breath and hold your breath for a few seconds, making sure the respirator collapses inward as you inhale.
3. Then exhale and hold your breath for a few seconds, making sure the respirator bulges outward as you exhale.

4. If air leaks, reposition the respirator, tighten the nose clip and reposition the elastic bands for a better fit. Follow steps 1-3 again until a tight seal has been achieved.
5. Respirators that have passed the fitting test are now safe to use.
6. Facial hair or other facial features may prevent the respirator from achieving a proper fit. If a tight seal is not achievable, please DO NOT use this product.
7. If you experience difficulty breathing, dizziness, or other symptoms, leave the contaminated area immediately.

## **WARNINGS**

1. The respirator does not eliminate the risk of contracting any disease or infection.
2. Improper use may lead to illness and even death.
3. Use this product immediately after package is opened.
4. Not recommended for more than 8 hours of use.
5. DO NOT sleep while wearing the respirator.
6. Avoid hand contact within the interior part of the respirator.
7. DO NOT use respirators if they expired.
8. For one-time use only. Dispose the respirator according to regulations.
9. DO NOT use if package is damaged.

## **Time use limitation**

If the respirator becomes damaged, soiled, or breathing becomes difficult, leave the contaminated area and replace and refit the new respirator.

## **Important notice**

1. To the extent permitted by law, A&Z shall not be liable for any loss or damage including any loss of business, loss of profits, or for any indirect, special, incidental or consequential loss or damage arising from reliance upon any information herein provided by A&Z. Nothing in this statement will be deemed exclude or restrict A&Z's liability for death or personal injury arising from its negligence.
2. Respirators are to be inspected prior to each use to assure there are no holes in the breathing zone.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

Manufactured by A&Z Pharmaceutical, Inc.  
350 Wireless Blvd., Hauppauge, NY 11788  
TEL: 1 (800) 480-1059

Made in U.S.A.