



## User Instructions for Respirator



- Before using, check the respirator for holes, tears or other damage to its surface, strap and nose cushioning to ensure that your respirator works well. Replace with a new respirator, if damaged.

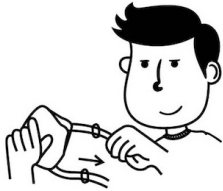

- Open the respirator and gently extend the elastic strap.
- Fold the nose clip into a  shape for a proper fit.



- Pull the lower strap over your head and place around the back of neck below your ears.


- Place the respirator against your face, with the bottom under your chin and nose cushioning across bridge of your nose.

- Pull the upper strap over your head and place high at the top back of your head.

- Lift the head strap a little for easier adjustment.

- Pull the strap outwards at either side of the buckle for tighter fit; push out the buckle for looser fit.

- Adjust the elastic strap for a proper and comfort tension.

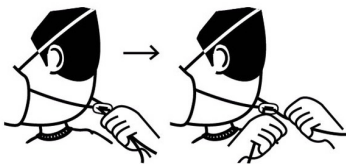

- Bend the nosepiece with both hands to fit snugly against your nose and face, and adjust respirator's position to ensure tight fit around face.



- To check the tight seal on your face, cover the middle surface of the respirator with one or both hands and exhale deeply. Do not push the respirator against your face. If you feel air blowing at your face or eyes, the respirator's position and strap need to be re-adjusted. Repeat the above steps until you do not feel air blowing at your face or eyes.