

User Instructions:

Doc id: *UIDBFFR05 v:2*

Filtering Facepiece Respirator : 203x series

*User Instructions for Dobu Mask Model
203V Filtering Facepiece Respirators*

Important information is provided in the User Instructions with each of the air filtration respirators 203x series that are NIOSH certified. **Failure to follow User Instructions for the air filtration being used may result in sickness or death.**



Warning

Misuse may result in sickness or death.

For proper use, see supervisor, or *User Instructions*, or call DOBU LIFE TECH in USA 3530 Wilshire Blvd St 1414, Los Angeles, 90010. +1 314 743 6243

Important

Before use, the user must read and understand the User Instructions. Keep these User Instructions for reference.

Use For

Respiratory protection from certain airborne contaminants according to NIOSH approvals, OSHA limitations,

User Instructions

- Failure to follow all instructions and limitations on the use of this respirator and/or failure to wear this respirator during all times of exposure can reduce respirator effectiveness and may result in sickness or death..
- Leave contaminated area immediately and contact supervisor if you smell or taste contaminants or if dizziness, irritation, or other distress occurs.
- Store respirator away from contaminated areas when not in use.
- Dispose of used product in accordance with applicable regulations

- Filtering facepieces are to be inspected prior to each use to assure there are no holes in the breathing zone other than the punctures around staples and no damage has occurred. enlarged holes resulting from ripped or torn filter material around staple/sonic welds punctures are considered damage. The User Instructions are also to provide information that the respirator has been tested and small punctures around the staples/sonic welds are normal and do not interfere with the respirator compliance with Part 84 approval requirements.

GENERAL SAFETY INFORMATION

Intended Use

The Dobu Mask 203x series filtering facepiece respirators are NIOSH approved and designed to help provide respiratory protection against certain airborne contaminants when used in accordance with all use instructions and limitations and applicable safety and health regulations.

Before occupational use of this product a written respiratory protection program must be implemented meeting all the requirements of OSHA 29 CFR 1910.134 such as training, medical evaluation, and fit testing and applicable OSHA substance specific standards. In Canada, CSA standard Z94.4 requirements must be met.

Use Limitations

- This respirator does not supply oxygen. Do not use in atmospheres containing less than 19.5% oxygen
- Do not use in atmospheres immediately dangerous to life or health
- Do not use when concentrations of contaminants are immediately dangerous to life and health, are unknown or when concentrations exceed 10 times the permissible exposure limit (PEL) or according to specific OSHA standards or applicable government regulations, whichever is lower
- Do not use when the U.S. Occupational Safety and Health Administration (OSHA) substance specific standards, such as those for, arsenic, cadmium or lead in the construction industry, specify other types of respiratory protection.
- Do not use with facial hair that can potentially prevent a good seal between the face and the face seal of the respirator
- Do not use when you have health problem that can hinder breathing.
- Do not use this product when your skin is known to have an allergic reaction, skin rash, or other irritant reaction to the respirator face piece

Time Use Limitation


If respirator becomes damaged, soiled or breathing becomes difficult, leave contaminated area immediately and replace the respirator

Storage Conditions and Shelf Life

Before use, store respirators in the original packaging and away from contaminated areas, dust, sunlight, extreme temperatures, excessive moisture and damaging chemicals. When stored in accordance with temperature and humidity conditions specified below, the product may be used until the "use by" date specified on packaging. Always inspect product and conduct a user seal check before use as specified in the User Instructions. If you cannot achieve a proper seal, do not use the respirator.



End of Shelf Life
Use respirators before the "use by" date specified on packaging

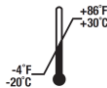

Storage Temperature Range
-20°C (-4°F) to +30°C (+86°F).



Storage Maximum Relative Humidity
<80% RH

Fitting Instructions

1. Hold the mask in your hand with the head straps hanging down.
2. Position the mask under your chin and over your nose. Slip the first strap over your head to the back of your neck.
3. Pull the upper strap into position at the back of your head, ensure the strap sits above your ears.
4. Adjust the nose clip to minimize leakage.
5. Perform a user seal check prior to each wearing. To check the respirator-to-face seal, place both hands completely over the respirator and inhale sharply. Be careful not to disturb the position of the respirator. A negative pressure should be felt inside the respirator. If any air leakage at the side of the mask, readjust the mask until it fits properly and ensure that the nose clip is positioned correctly on your nose.

If you cannot achieve a proper seal, DO NOT enter the contaminated area. See your supervisor.



Removal Instructions

See step 2 of Fitting Instructions and cup respirator in hand to maintain position on face. Pull bottom strap over head. Still holding respirator in position, pull top strap over head and remove respirator



Filter Configuration



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION

TC-	Protection	Respirator	Cautions and Limitations
		203V	
84A-DOB022a	N95	X	ABCJMNOP

1. Protection

N95 - Particulate Filter (95% filter efficiency level)
Effective against particulate aerosols free of oil;
time use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations"
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer"
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators."
- P - NIOSH does not evaluate respirators for use as surgical masks