


# FITTING INSTRUCTIONS



## Step 1

Hold respirator in your hand as shown, place the nosepiece at your fingertips and the headbands hanging freely below your hand. Inspect to make sure that there are no pin holes in the staple locations.

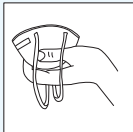

## Step 2

Place respirator under chin with the nosepiece facing up. Pull the top part of the headband to top of head.



## Step 3

Pull the bottom part of the headband over head and place below ears.



## Step 4

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.



## Step 5

Cup both hands over the mask, checking for tightness. Exhale sharply. Positive pressure should be felt in the mask. Re-adjust as necessary.

## Caution:

The respirator should be fit checked before wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved. Each mask is to be inspected prior to each use to assure there are no holes in the breathing zone other than the punctures around staples and no damage has occurred. Enlarged holes resulting from ripped or torn filter material around staple punctures are considered damaged.