



National Institute for Occupational Safety and Health



Instructions for Filtering Facepiece Respirator

Model #: GM-003-FFR-L / GM-005-FFR-S

Approval: TC-84A-PH04 / TC-84A-PH08

Assigned to: General Motors Company

P.O. BOX 33170

Detroit, MI 48232-5170

Protection: N95-Particulate Filter (95% filter efficiency level) Effective against particulate aerosols free of oil; time use restrictions may apply.

Precautions:


- A. Not for use in atmospheres containing less than 19.5 percent oxygen.
- B. Not for use in atmospheres immediately dangerous to life or health.
- C. Do not exceed maximum use concentrations established by regulatory standards.
- J. Failure to properly use and maintain this product could result in injury or death.
- L. Follow the manufacturer's User Instructions for changing cartridges, canister and/or filters.
- M. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O. Refer to User Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators for use as surgical masks.
- PH. The limited FFR approval will be effective until the declaration that circumstances exist justifying the authorization of the emergency use of personal respiratory protective devices during the COVID-19 outbreak is terminated under Section 564(b)(2) of the Act, or the EUA is revoked under Section 564(g) of the Act, or until such time that the approved configuration can no longer be maintained.

Prepare



Open respirator to extend with nosepiece in your hand.

Position

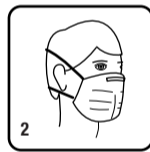

Position the respirator with nosepiece up. Pull straps over head with top at back and bottom around neck below ears.

Shape and Secure



Mold the nose area to the shape of your nose. Use both hands for effective respirator performance.

Seal Check



Place both hands completely over the respirator and inhale sharply. Be careful not to disturb the position of the respirator. A negative pressure should be felt inside the respirator. If air leaks around nose, readjust the nosepiece as described in step 3. If air leaks at the respirator edges, work the straps back along the sides of your head. If you CANNOT achieve the proper seal, DO NOT enter the contaminate area.

Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the requirements of OSHA 29 CFR 1910.134 such as training, fit testing, medical evaluation, and applicable OSHA substance specific standards.

Product Materials: Non-woven spunbond polypropylene fabric (outer layer) + Melt blown Polypropylene fabric (middle layers), Non-woven spunbond polypropylene fabric (inner layer).

Storage Conditions: Store in cool, ventilated and dry place.

For Concerns or Customer Feedback Email: facemask@gm.com