1730  
Respirator  
Disposable N95 Particulate Respirator and Surgical Mask

1730  
Respirator  
Disposable N95 Particulate Respirator and Surgical Mask

LD

|         | 1730 | 1740 |          |
|---------|------|------|----------|
| 844-160 | N95  |      |          |
| 844-161 | N95  | X    | ABCJKNOP |

1. Protection

N95 Particulate filter (95% filter efficiency level) effective against particulate aerosols free of oil, time use restrictions may apply.

2. Cautions and Limitations

A. Not for use in atmospheres containing less than 19.5% oxygen.  
 B. Not for use in atmospheres immediately dangerous to life or health.  
 C. Do not exceed maximum use concentrations established by regulatory standards.  
 J. Failure to properly use and maintain this product could result in injury or death.  
 M. All approved respirators shall be selected, fitted, used and maintained in accordance with NIOSH, OSHA, and other applicable regulations.  
 N. Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.  
 O. Refer to user instructions and/or maintenance manuals for information on use and maintenance of this respirator.  
 P. NIOSH does not evaluate respirators for use as surgical masks.

**IMPORTANT NOTICE TO PURCHASER**  
 Manufacturer warrants this product to be free from defects in material and workmanship. In the event that the product proves upon inspection to contain defects in material or workmanship, Manufacturer, at its option, to repair any portion showing defects or to refund to the Buyer the purchase price paid for this product. This warranty shall begin to run from the date of shipment of the product to the Buyer. This warranty shall be null and void if the product is used for any purpose other than that intended by the manufacturer. Except as provided herein, Manufacturer shall not be liable for any loss or damage, whether direct, consequential, or incidental, occasioned by the sale, use or misuse of this product, either then made or personal injury directly resulting from Manufacturer's negligence. Buyer shall determine the suitability of this product for its intended use and accept product shown, assumes all risks and liability in connection therewith.

LOUIS M. GERSON CO., INC.  
 Middleboro, MA 02348 U.S.A.  
 © 2014 All Rights Reserved. Made in U.S.A.  
 550488 rev 0  
 www.GersonCo.com

**Type N95 NIOSH Approved Respirator**

**WARNINGS**

- Read and understand before use.
- A properly selected respirator is essential to protect your health. Before using this or any product consult an Industrial Hygienist or Occupational Safety Professional to determine its suitability for your intended use.
  - This product will not protect against gases or vapors. Do not use in spray paint operations.
  - Do not use in sand blasting operations.
  - Do not use for protection against asbestos.
  - This product does not supply oxygen. Use this respirator in a well ventilated area containing sufficient oxygen to support life.
  - Do not use when concentrations of contaminants are unknown or immediately dangerous to life.
  - Leave work area immediately if: (A) breathing becomes difficult, or (B) dizziness or other signs of distress occur.
  - If user has any illness or disorders, consult a licensed medical doctor to determine suitability of use.
  - Facial hair and certain facial characteristics may prohibit effective use of this product.
  - Use only in accordance with instructions, labels and limitations relating to this product.
  - Do not alter or modify this product in any way.
  - Filtering facepieces are to be inspected prior to each use to assure there are no holes in the breathing zone and that no damage has occurred.

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Cat. No. 081730

1730  
Respirator  
Disposable N95 Particulate Respirator and Surgical Mask

FLUID RESISTANT  
 MEETS CDC GUIDELINES  
 FOR TB EXPOSURE CONTROL

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A surgical mask, when worn properly, may reduce potential exposure of the wearer to blood and bodily fluids, but does not eliminate the risk of contracting any disease or infection.

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

1. Hold the respirator in hand with nosepiece at your fingertips, allowing the headbands to hang freely below your hand.

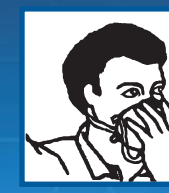
2. Press the respirator firmly against your face with the nose-piece on the bridge of your nose.

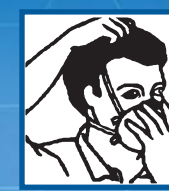3. Stretch and position the top band high on the back of your head. Stretch the bottom band over the head and position it below your ears.

4. Using both hands, mold the metal nosepiece to the shape of your nose.

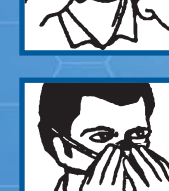5. To test fit cup both hands over the respirator and exhale vigorously. If air flows around your nose, tighten the respirator. If air leaks around the edges, reposition the straps for a better fit.

6. Remember, careful observance of these fitting instructions is an important step in safe respirator use.

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