


# N95 Respirators



**NIOSH, USA**  
**42CFR84**  
**Approved**

**TC-84A-6969**

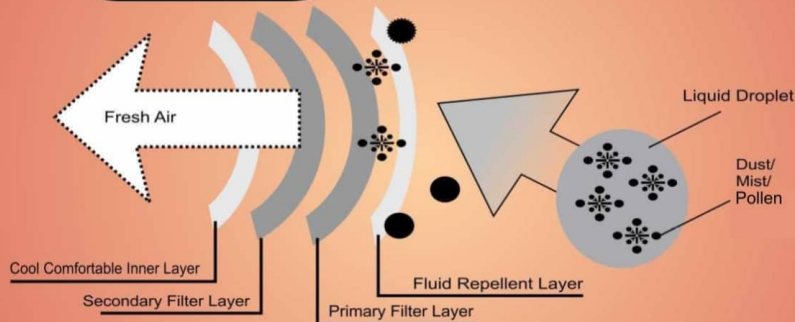

## User Instructions :

- 1) Place respirator under your chin with the nose piece up. pull the top head band over your head and above the ears.
- 2) Pull the bottom head band over your head and place around neck below the ears .
- 3) Push the nose piece inward with your finger tip to conform to the shape of your nose.

## User Seal Check Procedures

The individual who uses a tight-fitting respirator is to perform a user seal check to ensure that an adequate seal is achieved each time the respirator is put on. Either the positive and negative pressure checks listed in this appendix, or the respirator manufacturer's recommended user seal check method shall be used. User seal checks are not substitutes for qualitative or quantitative fit tests.

### I. Facepiece Positive and/or Negative Pressure Checks

A. Positive pressure check. Hold the mask with both hands & place it on your face and exhale gently into the facepiece. The face fit is considered satisfactory if a slight positive pressure can be built up inside the facepiece without any evidence of outward leakage of air at the seal.

Manufactured by :

**MAGNUM Health & Safety Pvt. Ltd.**

Formerly known as

**MAGNUM MEDICARE PVT. LTD.**

Survey No. 14, Hissa No.2, Unit No. 21,22,23 Supreme Estate, Village Devdal, Taluka Vasai,

Thane - 401208. Tel. : 022-28484282 • Customer Care No. : 08080850482 • e-mail : [info@magnumohs.com](mailto:info@magnumohs.com)

## N95 Respirators



B. Negative pressure check. Hold the mask with both hands & place it on your face, inhale gently so that the facepiece collapses slightly, and hold the breath for ten seconds. If the facepiece remains in its slightly collapsed condition and no inward leakage of air is detected, the tightness of the respirator is considered satisfactory.

### II. Manufacturer's Recommended User Seal Check Procedures

The respirator manufacturer's recommended procedures for performing a user seal check may be used instead of the positive and/or negative pressure check procedures provided that the employer demonstrates that the manufacturer's procedures are equally effective.

## Warnings :

- 1) A properly selected respirator is essential to protect your health. Before using this or any product, consult an Industrial Hygienist or Occupational Professional to determine its suitability for your intended use.
- 2) This product is for dry aerosols only. It will not protect against gases or vapours or in spray paint operations.
- 3) Do not use in sand blasting operations.
- 4) Do not use for protection against asbestos.
- 5) This product does not supply oxygen. Use only in well ventilated areas containing sufficient oxygen to support life.
- 6) Do not use when concentrations of contaminants are unknown or immediately dangerous to life, or when concentrations exceed 10X PEL.
- 7) Leave work area immediately if: (a) breathing becomes difficult; (b) dizziness or other signs of distress occur.
- 8) If user has any illness or disorders, consult a licensed medical doctor to determine suitability of use.
- 9) Facial hair and certain facial characteristics may prohibit effective use of the product.
- 10) Use only in accordance with instructions relating to this product.
- 11) Discard and replace respirator if it becomes damaged or breathing resistance becomes excessive.
- 12) Do not alter or modify this product in any way.
- 13) Follow the manufacturer's instructions for fitting the respirator to the face.

**NOTE : Before issuing any respirator to be worn, a face fit test must be performed.**



**MAGNUM HEALTH & SAFETY PVT. LTD.**



Correspondence: 602 Mamta Heights, Devidas Lane, Borivali(West) Mumbai-400 103, Maharashtra, India

Factory: Survey No.14, Hissa No.2, Unit No.21,22,23 Supreme Estate,Village Devdal, Taluka Vasai,Thane

Tel:022-65948415 Tel(fax):022-28484282 e-mail : [info@magnumohs.com](mailto:info@magnumohs.com) or [magnum.disp@gmail.com](mailto:magnum.disp@gmail.com)

These respirators are approved only in the following configurations:

TC-	Protection	Respirator	Cautions and limitations
84A-6969	N95	N95001	ABCJMNOP
		X	

### 1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols  
free of oil; time use restrictions may apply

### 2.Cautions and Limitations

A-Not for use in atmospheres containing less than 19.5 percent oxygen

B-Not for use in atmospheres immediately dangerous to life or health.

C-Do not exceed maximum use concentrations established by regulatory standards.

J-Failure to properly use and maintain this product could result in injury or death.

M-All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.

N-Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.

O-Refer to User's Instructions and / or maintenance manuals for information on use and maintenance of these respirators.

P-NIOSH does not evaluate respirators for use as surgical masks.

"Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United State employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing."

Manufactured by :

**MAGNUM Health & Safety Pvt. Ltd.**

Formerly known as

**MAGNUM MEDICARE PVT. LTD.**

Survey No. 14, Hissa No.2, Unit No. 21,22,23 Supreme Estate, Village Devdal, Taluka Vasai,

Thane - 401208. Tel. : 022-28484282 • Customer Care No. : 08080850482 • e-mail : [info@magnumohs.com](mailto:info@magnumohs.com)