

## FITTING INSTRUCTIONS:

1



### Step 1

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

2



### Step 2

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.


3



### Step 3

Pull the bottom headband overhead and place around neck below ears.


4

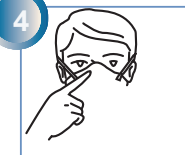

### Step 4

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.

5

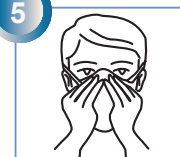

### Step 5

**Cup both hands over the mask, checking for tightness, Inhale sharply, negative pressure should be felt in the mask, Re-adjust as necessary.**

### Caution:

The respirator should be fit checked for each wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved.