



Form-N-Fit™ N95 Respirator Mask

Model Numbers:
110.118 / 119.128
123.138 / 107.128
124.124




Limited Warranty

This product is warranted to be free of defects in material and workmanship for a period of six (6) months from the date of original purchase. ThermoPore will, at its option, repair or replace without charge, except for transportation costs, parts that fail under normal use and service when operated and maintained in accordance with our instructions. This warranty does not cover damage caused by misuse, abuse, or unauthorized modifications. Proof of purchase is required to verify eligibility.

Use For: Solids such as those from processing materials, coal, iron, ore, flour, and certain other substances. Liquids or non-oil based dissolved particles from sprays that do not also emit oil aerosols or vapors.

Do Not Use: Gases and vapors, including those present in paint spraying operations oil aerosols, asbestos, or sandblasting.

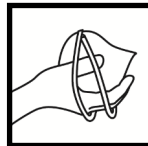

Step 1
Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.



Step 2
Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



Step 3
Stretch and position the top band high on the back of your head. Stretch the bottom band over the head and position below your ears.



Step 4
Using both hands, mold the nose piece to the shape of your nose.



Step 5
To test fit: Cup both hands over the respirator being careful not to disturb position, and inhale aggressively. If air leaks around edges, reposition the straps of adjust tension for better fit. Follow these instructions during each use to achieve proper fit.

Warnings

Wearer must read fitting and use limitations before use.

This respirator does not protect against the risk of contracting disease or infection.

Individuals highly sensitive to natural rubber latex may have an allergic reaction.

Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the requirements of OSHA 29 CFR 1910.134 such as training, fit testing, medical evaluation, and applicable OSHA substance specific standards.

If you CANNOT achieve a proper seal, DO NOT enter the contaminated area.

Service Time Limitation

Do not wear the same respirator for more than 8 hours, continuous or intermittent, in dirty workplace that could result in high filter loading.

Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted.

The maximum service time of this respirator is seven (7) consecutive calendar days (including day of non-use), beginning from the first day of use.



ThermoPore Materials Corp.
27 First Avenue
Suite #350
Newnan, GA 30263
678.621.0200



TC-84A-	Protection	Respirator	Cautions and Limitation
PH####	N95	110.118	A B C J M N O P P H
PH####	N95	119.128	A B C J M N O P P H
PH####	N95	123.138	A B C J M N O P P H
PH####	N95	107.128	A B C J M N O P P H
PH####	N95	124.124	A B C J M N O P P H

PROTECTION

N95 - Particle Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

CAUTIONS and LIMITATIONS

- A Not for use in atmospheres containing less than 19.5 percent oxygen.
- B Not for use in atmospheres immediately dangerous to life or health.
- C Do not exceed maximum use concentrations established by regulatory standards.
- J Failure to properly use and maintain this product could result in injury or death.
- M All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O Refer to User Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P NIOSH does not evaluate respirators for use as surgical masks.
- PH The limited FFR approval will be effective until the declaration that circumstances exist justifying the authorization of the emergency use of personal respiratory protective devices during the COVID-19 outbreak is terminated under Section 564(b)(2) of the Act, or the EUA is revoked under Section 564(g) of the Act, or until such time that the approved configuration can no longer be maintained.