


FITTING INSTRUCTION -

使用方法

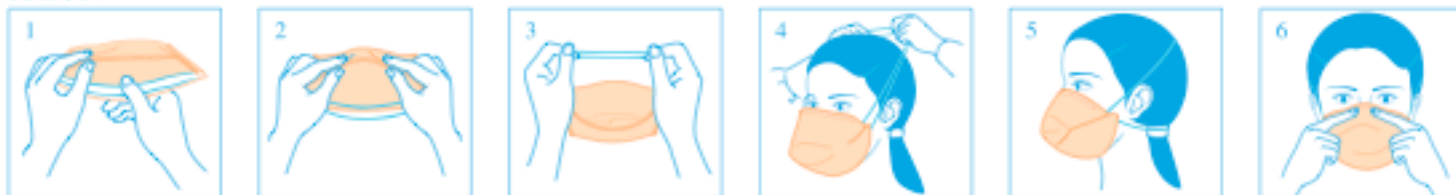

1. Separate mask to open.
2. Bend the nosepiece to conform mask to your nose.
3. Hold the center of headbands separately.
4. Hang the lower part on your chin and bring both bands over your head.
5. Place the upper band at your crown and the lower band at your neck.
6. Form the nosepiece tightly to make no chink.
※ Using both hands,over the whole of the respirator to check the fit-breathing out to make sure that no breath escapes from around the nose, cheeks or edges of the respirator.

1. ノーズクランプ側を上にし、マスクを上下に開きます。
2. ノーズクランプを鼻に合うように軽く曲げます。
3. 2本のゴムを別々に持ちます。
4. マスク下側をアゴに掛け、2本のゴムの中央を持ち、頭の後へ持っていきます。
5. マスク上側ゴムを頭頂部に、下側ゴムを首に固定します。
6. すき間ができないようにしっかりとノーズクランプを顔に合わせます。
※ 最後にマスクが密着しているか確認します。
先ず、両手でマスクを覆い、息を吐き出し、空気のもれをチェックして密着の良い状態にマスクを合わせます。