

Venus CN95

Particulate Respirator N95

USE INSTRUCTIONS :


1. Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the whole duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury or death.
2. Before occupational use, the wearer must be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards.
3. OSHA standards including 29CFR1910.134 (e) (5) and CSA standard Z 94.4-93 require that the wearer be fit-tested
4. Discard the respirator and replace with a new one if:
 - Excessive clogging of the respirator causes breathing difficulty
 - The respirator becomes damaged.Discard the respirator and replace with a new one after every use.
5. Leave the contaminated area if dizziness, irritation or other distress occurs.
6. Do not use with beards or other facial hair or other conditions that prevent a good fit between the face and the fitting surface of the respirator.

STORAGE INSTRUCTIONS :

The filter respirator until use shall be stored in the sealed single pack to retain its properties. For transport such singles packs shall be suitably packed in outer cartons to protect from climatic hazards & mechanical shocks.


SHELF LIFE :

60 months from the date of manufacture (if stored between -5°C - $+50^{\circ}\text{C}$ & Humidity not over 80%)

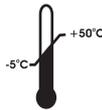


DONNING AND FITTING INSTRUCTIONS:

Follow these instructions each time respirator is worn.

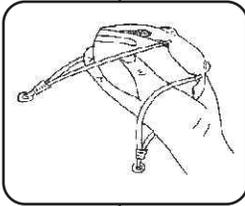


- Remove respirator from package and “pre-form” the aluminum nose bridge adjuster.

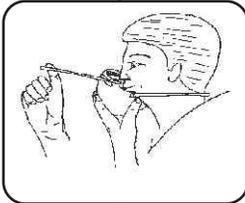


- Hold the respirator in the hand, allowing the head-straps to hang loosely in front of the mask.

- Position the lower head strap around the neck below the ears. Hold the respirator against the face and place the top head-strap around the crown of the head.

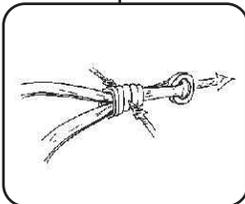

- Adjust head straps by grasping the locking slide and pull adjuster ring until respirator fits comfortably.

- Fit the nose area to the shape of the face by forming the nose piece with the fingertips. Start from the top of the nose piece working downward while pressing against the nose. Always use both hands.

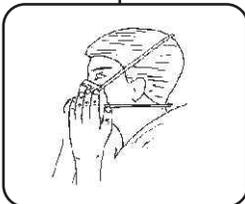

FIT-CHECK THE FACE-SEAL AS FOLLOWS:

- Place both hands over the respirator - do not disturb the position of the respirator.
- Exhale sharply to detect leakage.

WARNING:

If you cannot achieve satisfactory fit, do NOT enter the contaminated area.
See your Supervisor



Venus Safety & Health Pvt. Ltd.
Plot No. L-17, MIDC, Taloja, Tondra
Navi-Mumbai : 410 208. Maharashtra, India
E-mail: info@venusohs.com Web: www.venusohs.com
Telephone no: +91-22-27401663 Fax no: +91-22-27411652

THESE RESPIRATORS ARE APPROVED ONLY IN THE FOLLOWING CONFIGURATIONS:

TC-	Protection ¹	Respirators	Caution and Limitations ²
		CN95	
84A-5492	N95	X	ABCJMNOP

1. Protection:

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil, time use restrictions may apply

2. Cautions and Limitations:

- A - Not for use in atmosphere containing less than 19.5 % oxygen.
- B - Not for use in atmosphere immediately dangerous to life or health.
- C - Do not exceed maximum use concentration established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.