

VENUS CN95 ME -

Particulate respirator N95 -

USE INSTRUCTIONS :

1. Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the whole duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury or death.
2. Before occupational use, the wearer must be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards.
3. OSHA standards including 29CFR1910.134 (e) (5) and CSA standard Z 94.4-93 require that the wearer be fit-tested.
4. Discard the respirator and replace with a new one if:
 - Excessive clogging of the respirator causes breathing difficulty
 - The respirator becomes damaged. -Discard the respirator and replace with a new one after every use. -
5. Leave the contaminated area if dizziness, irritation or other distress occurs.

USE LIMITATIONS:

1. Do NOT use the respirator or enter or stay in a contaminated area under the following circumstances:
 - Atmosphere contains less than 19.5% oxygen.
 - If you smell or taste a contaminant
 - For protection against gases or vapors
 - Contaminants or their concentration are unknown or immediately dangerous to life or health.
 - Construction of contaminants exceed maximum use constructions in applicable OSHA standards or applicable government regulations or 10 times the PEL (Permissible Exposure Limit), whichever is lower.
 - For sandblasting, paint-spray operations, asbestos.
2. Do NOT modify or misuse the respirator.
3. Do NOT use the respirator with beard or other facial hairs that interferes with direct contact between the face and the edge of the respirator, or any other conditions that may prevent a good face-seal.

“A” Series respirators should be replaced when the wearer senses an increase in breathing resistance due to the filters collecting particulate matter. Replace mask frequently.

PROTECTION :

A-series filters: effective against particulate aerosols free of oil.

Do NOT use

- For protection against hazardous concentration of gases or vapor
- For sand-blasting, paint-spray operations, asbestos

FILTERING FACEPIECE APPROVAL LABEL



Mfg. By :
Venus Safety & Health Pvt. Ltd.
Plot No. L-17, MIDC, Taloja, Tondra
Navi Mumbai - 410 208, Maharashtra, India
E-mail : info@venusohs.com Web. : www.venusohs.com
Telephone no: +91-22-27401663 Fax no: +91-22-27411652



THESE RESPIRATORS ARE APPROVED ONLY IN THE FOLLOWING CONFIGURATIONS:

TC-	Protection ¹	Respirators		Cautions & Limitations ²
		CN95	CN95 ME	
84A-5492	N95	X		ABCJMNO P
84A-5634	N95		X	ABCJMNO P

1. Protection:


N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil, time use restrictions may apply

2. Cautions and Limitations:


- A - Not for use in atmospheres containing less than 19.5 % oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.

DONNING AND FITTING INSTRUCTIONS:


Follow these instructions each time respirator is worn.

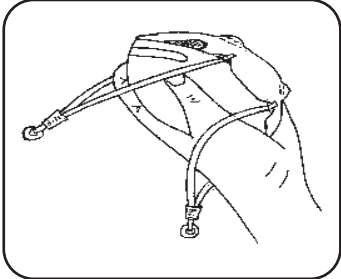

- Remove respirator from package and “pre-form” the aluminum nose bridge adjuster.

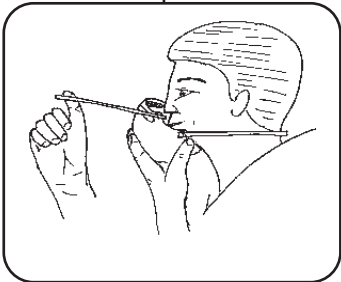


- Hold the respirator in the hand, allowing the head-straps to hang loosely in front of the mask.
- Position the lower head strap around the neck below the ears. Hold the respirator against the face and place the top head-strap around the crown of the head.

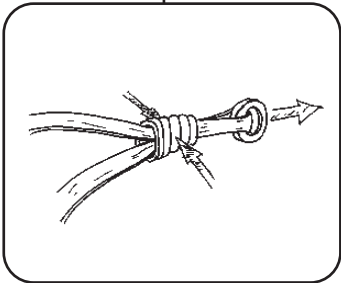

- Adjust head straps by grasping the locking slide and pull adjuster ring until respirator fits comfortably.

- Fit the nose area to the shape of the face by forming the nose piece with the fingertips. Start from the top of the nose piece working downward while pressing against the nose. Always use both hands.

FIT-CHECK THE FACE-SEAL AS FOLLOWS:

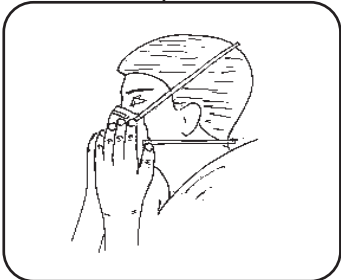

- Place both hands over the respirator - do not disturb the position of the respirator.
- Exhale sharply to detect leakage for unvalved respirator. - (Condition does not apply for valved respirator)
- Inhale sharply to detect leakage for valved respirator. - (Condition does not apply for unvalved respirator) -

WARNING:

If you cannot achieve satisfactory fit, do NOT enter the contaminated area. See your Supervisor.

This respirator helps protect against certain contaminants. It may not eliminate the risk of contracting disease or infection. Misuse may result in sickness or death.

FOR MORE INFORMATION or assistance on this product, contact :

VENUS (Technical Service Dept)

+91-22-27692646 e-mail : info@venusohs.com, www.venusohs.com. -