


Donning & Fitting Instruction



1. Adjust the elastic straps for suitable size before placing the respirator on the face.

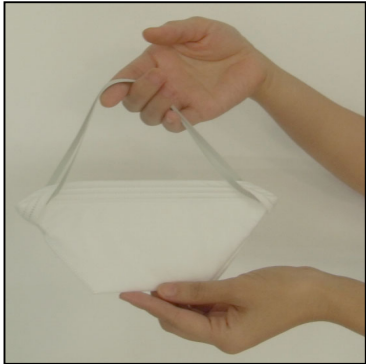

2. Cup the respirator in your hand, with nose piece towards your fingers.



3. Position the respirator at mouth and nose. Pull the elastic strap near your hand over your head resting it at the back of your head. Pull the strap near your finger over your head and position it at the top of your head and adjust the elastic straps.



4. Shape the nose piece suitably with your forefingers and middle fingers in order to keep the respirator to cling your face.



5. Use two hands perform a user seal check prior to each wearing. To check fit, place both hands completely over the respirator and exhale. If air leaks around nose, readjust the nosepiece as described in step 4. If air leaks at the respirator edges, work the straps back along the sides of your head.

Removal Instructions

Cup respirator in hand to maintain position of face and pull bottom strap over head. Still holding the respirator in position, pull top strap over head and remove respirator.

Dispose of used product in accordance with applicable regulations.

Rev. A6. June. 1, 2024