

FILTERING FACEPIECE APPROVAL LABEL



Latoplast Ltd.
1670 Finfar Ct, Unit 1
Mississauga, ON L5J4K1
1-800-268-5340



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:



TC-	Protection ¹	Respirator	Cautions and Limitations ²
		027-44150	
84A-7719	N95	X	ABCJMNOPS

1. Protection

N95- Particulate Filter (95% filter efficiency level)
effective against particulate aerosols free of oil;
time use restrictions may apply

2. Cautions and Limitations

- A- Not for use in atmospheres containing less than 19.5% oxygen.
- B- Not for use in atmospheres immediately dangerous to life or health.
- C- Do not exceed maximum use concentrations established by regulatory standards.
- J- Failure to properly use and maintain this product could result in injury or death.
- M- All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N- Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O- Refer to User Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P- NIOSH does not evaluate respirators for use as surgical masks.
- S- Special or critical User Instructions and/or specific use limitations apply. Refer to User Instructions before donning.

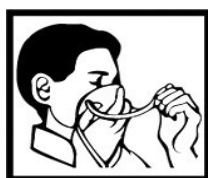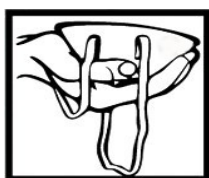

Instructions for the 027-44150 N95 Maintenance free respirator Fitting Instructions

The respirator user must be fit tested according to OSHA regulation 29 CFR 1910.134.

1. Hold the respirator in one hand with the nose area at your fingertips, allowing the headbands to hang below your hand.
2. Press the respirator against your face with the nose foam on the bridge of your nose.
3. Place the top band high on the back of your head. Move the bottom band over your head and position it below your ears.
4. Test the fit. Cup both hands over the respirator and exhale vigorously. If any leakage is detected, readjust the position of the respirator according to steps #2, #3, and #4 till get a better fit.

If you CANNOT achieve a proper fit, DO NOT enter the contaminated area.
See your supervisor for proper fitting respiratory protection.

APPLICATIONS

USE AGAINST:

NIOSH approved for protection against particulate aerosols free of oils.

Latoplast Ltd. recommends this respirator for use against particulates such as flour, iron coal and minerals.

DO NOT USE AGAINST:

1. Vapors, paint spray, oil-based aerosols, asbestos and sand blasting.
2. Any contaminant whose concentration is unknown, or immediately dangerous to life.
3. Concentrations of Dust and Mists which exceed the maximum use concentration or 10 times contaminant's permissible exposure limit established by OSHA.

WARNING:

1. These respirators do not supply oxygen. Use only in well ventilated areas containing sufficient oxygen to support life.
2. Leave contaminated area immediately if the respirator becomes damaged and proper fit cannot be obtained or breathing becomes difficult. Dispose and change the respirator.

IMPORTANT LIMITATIONS:

1. All instructions and warnings have to be followed. Failure to do so will reduce the efficiency of respirator and cause sickness or death.
2. Training on the proper use of respirator must be provided by the employer in accordance with the safety and health standard.
3. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
4. Do not alter or modify the respirator.
5. Not suitable for persons with beard or other facial hair which may reduce the efficiency of respirator.
6. The use and reuse of this respirator is subject to consideration of hygiene damage and increased breathing resistance. Service time should not be extended beyond 8 total hours of use.
7. Store the respirators in cool place. Avoid direct exposure of sunlight during storage as the filter media in this respirator can be degraded by extended exposure to direct sunlight.
8. S-Special or critical User Instructions: The model 027-44150 filter facepiece respirator has been manufactured by JINFUYU INDUSTRIAL CO., LTD. for Latoplast Ltd. under TC-84A-7719.

N95: 95% filtration efficiency against solid and liquid aerosols that do not contain oil.

ATTENTION: READ APPLICATIONS, LIMITATIONS AND FITTING INSTRUCTIONS ON BOX BEFORE RESPIRATOR IS WORN.

CONTENTS: 20 RESPIRATORS

