



**FILTERING FACEPIECE
APPROVAL LABEL**



**TUNGSAFE CORPORATION
(BRAND NAME: BLUE EAGLE)**

**139 HO-PING 2nd Road, Kaohsiung 80611, Taiwan
Tel: +886-7-722-6666**

THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	Protection¹	Respirator	Cautions and Limitations²
		F750C	
84A-7836	N95	X	ABCJMNOPS

1. Protection

N95 - Particulate Filter (95% filter efficiency level)
Effective against particulate aerosols free of oil;
time use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.
- S - Special or critical User's Instructions and/or specific use limitations apply.
Refer to User's Instructions before donning.

FIDO MASKS Co., LTD.
USER INSTRUCTIONS
F750C N95 Particulate Respirator

User's Instructions



WARNING TO USER

Before use, the user must read, understand and follow all User's Instructions, warnings and use limitations of this respirator. This respirator provides limited protection and must be properly selected and fit tested by a competent person to provide protection from known contaminants according to OSHA 29 CFR 1910.134 or applicable safety and health standard. The respirator's face seal must be checked each time before entering a contaminated area and the respirator must be worn during all times of exposure to contaminants. **Failure to follow all of these instructions, warnings and limitations, contained on this package may reduce the effectiveness of this respirator and result in serious sickness or DEATH!!** This respirator offers nuisance level relief from (type of odor (such as organic vapors)) that are below the Permissible Exposure Limit (PEL). Nuisance level refers to concentrations not exceeding the OSHA PEL or other government occupational exposure limits, whichever is lower.

The F750C N95 respirator with a filtration efficiency of at least 95% against solid or liquid particles that do not contain oil.

DO NOT USE FOR:

The F750C N95 respirator should not be used for gases, vapors, oil aerosols or where concentrations of contaminants exceed OSHA or applicable safety and health standard's maximum use concentrations, are unknown or immediately dangerous to life and health.

USE INSTRUCTIONS:

1. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
2. Prior to each use, carefully inspect the entire respirator for correct assembly and for tears or damage including the straps, filter and exhalation valve (if equipped with one). ***Discard the entire respirator if any part does not function properly*** and get a new respirator.
3. If the respirator becomes damaged, distorted, soiled, if breathing becomes difficult or a proper respirator to face seal cannot be maintained, ***leave the contaminated area immediately*** and get a new respirator.

4. If dizziness, irritation, nausea or other distress occurs, ***leave the contaminated area immediately*** and contact your supervisor or medical professional.
5. Do not alter abuse or misuse this respirator.
6. Properly dispose of used respirators in accordance with applicable regulations.

USE LIMITATIONS:

1. This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.
2. This respirator is not to be used when the user has beards, facial hair, scars or other conditions that prevent a good seal between the face and the respirator.
3. This respirator is designed for occupational or professional use by adults who are properly trained in their use and limitations. Respirators must not be used by children.
4. Particulates and other contaminants which can be dangerous to your health include those that you cannot see, taste or smell. Respirators can help to protect your lungs against certain airborne contaminants, but cannot prevent entry through other routes such as the skin. Additional types of personal protective equipment may be required.

STORAGE:

1. Store unused respirators in their sealed, original packaging in a clean, dry, non-contaminated area that is protected from extreme temperatures, moisture, sunlight, dust and chemicals.
2. Never attempt to clean and reuse this respirator.
3. ***Used respirators should not be stored*** and should be disposed of properly.


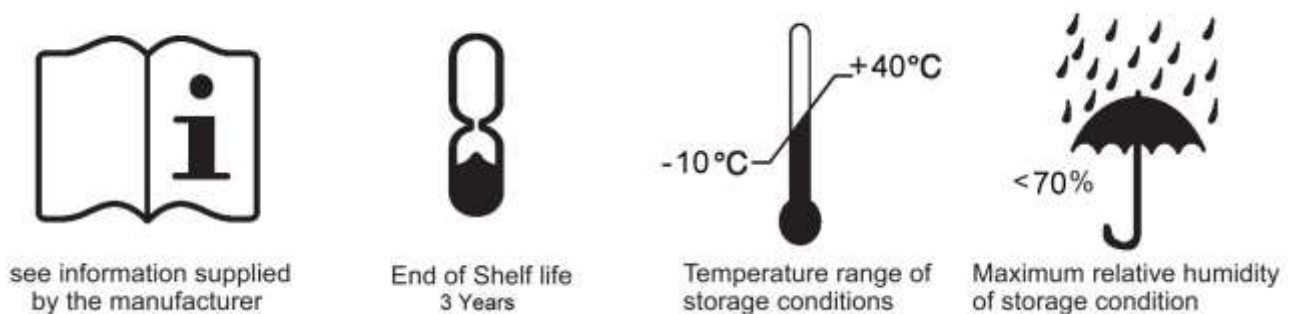


Fig.3 Box printing side 1

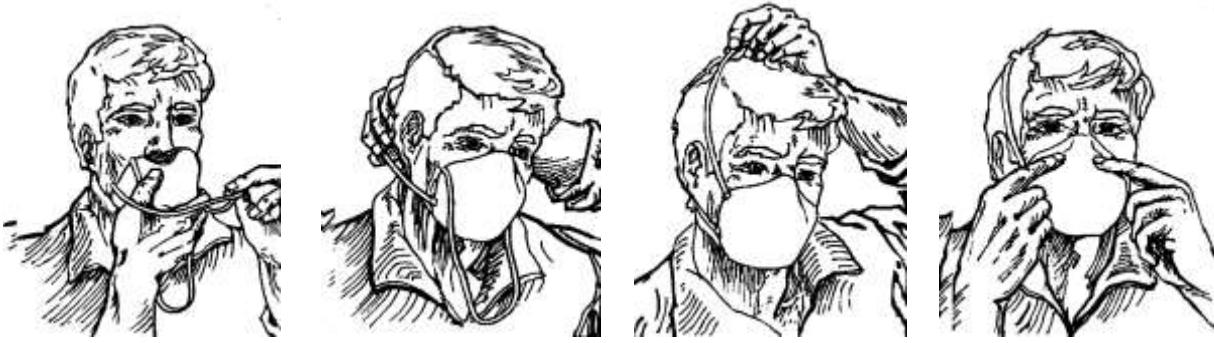

FITTING INSTRUCTIONS:

Before fitting the respirator, *carefully inspect the entire respirator* to ensure that the respirator is not damaged and is functioning correctly. *Discard the entire respirator if any part does not function properly.*

1. Hold the respirator in your hand with the nosepiece at your fingertips and allowing the straps to fall below the hand.
2. Place respirator under the chin with the nosepiece (narrow end) up. Pull the continuous bottom strap over your head and position it around the neck below the ears.
3. Raise the continuous top strap over your head and place it high to the back of the head.
4. Using both hands, push the nosepiece inwards with your fingertips to mold and shape the nose area to the face.
5. A seal check **MUST** be performed each time the user enters a contaminated area. Lightly cup both hands over the front of the respirator and **INHALE SHARPLY**. A negative pressure should be felt inside of the respirator. If any air flow is detected around the edges of the respirator, readjust the straps. Repeat seal check until a proper seal is achieved; otherwise, see your supervisor. ***NEVER enter a contaminated area without a proper seal.***

Note: Pinching the nosepiece using one hand may result in improper fit and less effective respirator performance. Use two hands.

Special Instructions: The model F750C filtering facepiece respirator has been manufactured by Approval Holder Fido Masks CO., Ltd. for TUNGSAFE CORPORATION under TC-84A-7836.

Lot No. XXXXXX

(the Lot No. XXXXXX will be as same as abbreviated Label indicated, the "XXXXXX" should be printed as purchase order number)

Contents: 20 pcs.

Fig. 4 Box printing side 3