

# Model FT-N142 Particulate Respirator N95

*User Instructions*  
*(Keep insert for reference)*

Rev. A  
08/15/2021




## Warning

This respirator helps protect against certain particles. **Misuse may result in sickness or death.** For proper use, see supervisor, *User Instructions*, or call +86-512-65586028 / 65586078

Company Name: Suzhou Fangtian Industries Co., Ltd.  
Address: Xiehe Building, No.158 Jinmen Road, Suzhou City, Jiangsu, P. R. China  
Fax: +86-512-65586038  
E-mail: [sales@chinamasks.com](mailto:sales@chinamasks.com)  
<http://www.chinamasks.com/>

### USE LIMITATIONS:

1. Do **NOT** use the respirator or enter or stay in a contaminated area under the following circumstances:
  - Atmosphere contains less than 19.5% oxygen
  - Atmosphere contains oil aerosols if using N protection class respirator
  - For protection against gases or vapors
  - Contaminants or their concentrations are unknown or immediately dangerous to life or health
  - Concentrations or contaminants exceed maximum use concentrations in applicable OSHA standards or applicable government regulations or 10 times the PEL (Permissible Exposure Limit), whichever is lower
  - For sandblasting, paint-spray operations, asbestos
2. Do **NOT** modify or misuse the respirator.
3. Do **NOT** use the respirator with beards or other facial hair that interferes with direct contact between the face and the edge of the respirator, or any other conditions that may prevent a good face-seal.

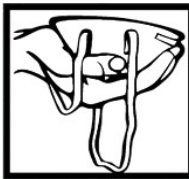

### WARNING:

This respirator helps protect against certain contaminants. It may not eliminate the risk of contracting disease or infection. **Misuse may result in sickness or death.**

Materials which may come into contact with the wearer's skin could cause allergic reactions to susceptible individuals.

FOR MORE INFORMATION or assistance please contact with Manufacturer.

### FITTING INSTRUCTIONS

Follow these instructions each time respirator is worn.

1. Hold the respirator in one hand with the nosepiece at your fingertips, allowing the headbands to hang below your hand.
2. Press the respirator against your face with the nosepiece on the bridge of your nose.
3. Place the top band high on the back of your head. Move the bottom band over your head and position it below your ears.
4. Using both hands, mould the nosepiece to the shape of your nose.
5. Test the fit. Cup both hands over the respirator and exhale vigorously. If air flows around your nose, tighten the nosepiece. If air leaks around the edges, reposition the bands for better fit.
6. Fit-check the face-seal as follows:
  - Place both hands over the respirator without disturbing its position.
  - If the respirator is not fitted with an exhalation valve, exhale sharply.
  - If the respirator is fitted with an exhalation valve, inhale sharply (a negative pressure should be felt inside the respirator).
  - If you detect air-leaks, readjust the head straps and/or the position.

**WARNING:** If you cannot achieve satisfactory fit, do NOT enter the contaminated area. See your Supervisor.

### USE INSTRUCTIONS:

1. Failure to follow all instructions and limitations concerning the use of the respirator and/or failure to wear the respirator for the whole duration of exposure to contaminants can seriously reduce the performance characteristics of the respirator and lead to illness, injury or death.
2. Before use, the wearer must first be trained by the employer in the correct use of the respirator in accordance with applicable safety and health standards.
  3. OSHA standards including 29 CFR 1910.134(e) (5) and CSA standard Z94.4-93 require that the wearer be fit-tested. Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
4. Discard the respirator and replace with a new one if:
  - Excessive clogging of the respirator causes breathing difficulty.
  - The respirator becomes damaged.
5. Leave the contaminated area if dizziness, irritation or other distress occurs.





**Suzhou Fangtian Industries Co., Ltd.**

Office: Xiehe Building, No.158 Jimmen Road, Suzhou City, Jiangsu, P. R. China

Factory: Xingfu South Rd, Medical Products Industrial Park, Yandu District, Yancheng, City, Jiangsu Province

TEL :+86-512-65586028 FAX : +86-512-65586038

THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
		FT-N142	
<b>84A-9425</b>	N95	X	ABCJMNO P

**1. Protection**

N95 - Particulate Filter (95% filter efficiency level)  
Effective against particulate aerosols free of oil;  
time use restrictions may apply

**2. Cautions and limitations**

A-Not for use in atmospheres containing less than 19,5 percent oxygen.

B-Not for use in atmospheres immediately dangerous to life or health.

C-Do not exceed maximum use concentrations established by regulatory standards.

J-Failure to properly use and maintain this product could result in injury or death.

M-All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.

N-Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.

O-Refer to User's Instructions and/or maintenance manuals for information on use and maintenance of these respirators.

P-NIOSH does not evaluate respirators for use as surgical masks.

**NIOSH APPROVED: N95**

95% efficient against particulate aerosols free of oil. This respirator helps protect against certain particulate contaminants but does not eliminate exposure to or the risk of contacting disease or infection. Misuse may result in sickness or death.