




# AP0028+ SURGICAL N95 RESPIRATOR FITTING INSTRUCTION



**Step 1:**  
Fold the respirator  
from right to left.

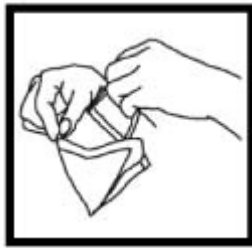

**Step 2:**  
Pull the both straps  
out from the folded  
respirator.



**Step 3:**  
Position the respirator  
under the chin with  
noseclip facing upwards.



**Step 4:**  
Pull the straps over  
your head and rest  
them on your neck.



**Step 5:**  
Pull the one of the  
strap over your head  
and position it.



**Step 6:**  
Mould the noseclip to  
the shape of the nose.

- Change respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted, or a proper face fit cannot be maintained. Remember, careful observance of these fitting instructions is an important step in safe respirator use.
- To check the correct fit, cup with both hands the respirator and exhale powerfully. If any airflow indicates leakage around the edges of the respirator, stretch the headbands again by pulling them backwards. Then correct the position of the respirator (see picture) and check again.

**WARNING!** It is not allowed to use the respirator without proper fit. If you can't achieve a proper seal, do not use the respirator.