


FITTING INSTRUCTION

STEP 1 : Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

STEP 2 : Open and place respirator on chin. Stretch upper strap and position high on back of head. Stretch bottom strap over the head and position below ears. Using both hands, mold the metal nosepiece to the shape of your nose bridge.

STEP 3 : To fit check a) cup both hands over respirator be careful not to disturb position, b) inhale vigorously, c) exhale vigorously. If air leaks around the edges during inhalation or exhalation, re-position straps or adjust straps tension for better fit. Repeat (a), (b) and (c). If there are still noticeable leakage then this respirator may not be suitable for you and we advise you to seek another style from our range.

Careful observance of these fitting instructions is an important step in safe respirator use.

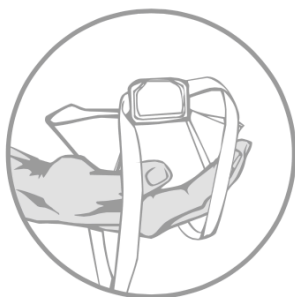

The respirator should be fit checked prior to each wearing.