


Donning Instructions:

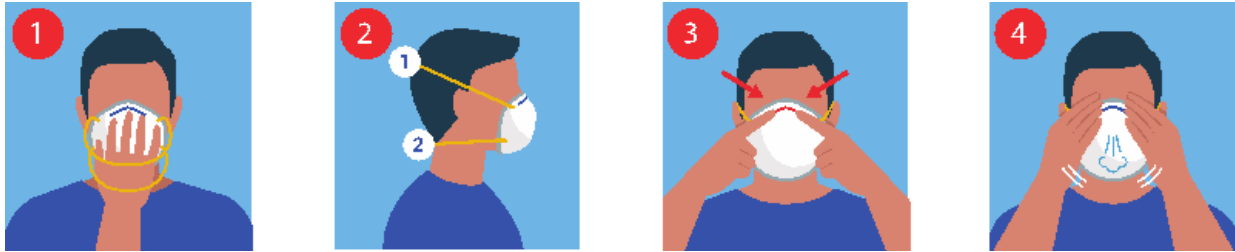

1. Open it into a cup shape and cup the respirator in your hand, with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.
2. Position the respirator under your chin with the nosepiece up. Pull the top strap over your head resting it high at the top back of your head. Pull the bottom strap over your head and position it around the neck below the ears. Make certain hair, facial hair, jewelry, and clothing are not between your face and the respirator as they will interfere with fit.
3. Place your fingertips from both hands at the top of the metal nosepiece. Using two hands, mold the nose area to the shape of your nose by pushing inward while moving your fingertips down both sides of the nosepiece. **Note:** Pinching the nosepiece using one hand may result in improper fit and less effective respirator performance (Use two hands).
4. Perform a User Seal Check. To check the respirator-to-face seal, place both hands completely over the respirator and exhale sharply. Be careful not to disturb the position of the respirator. If air leaks around nose, readjust the nosepiece as described in step 3. If air leaks around the respirator edges, adjust position of straps, and make certain respirator edges fit snugly against the face.

If you CANNOT achieve a proper seal, DO NOT enter any contaminated areas. See your supervisor.

Removal Instructions:

Without touching the respirator, slowly lift the bottom strap from around your neck up and over your head. Then lift off the top strap. Dispose of used product in accordance with applicable regulations. **WASH YOUR HANDS AFTER USE!**