

# Surgical/N95 Filtering Facepiece Respirator (RES-12-A / RES-12-B / RES-12-C)

## User Instructions

Document ID: UI-02

Revision B

### Use for:

Reduction of exposure to particles such as dust, pollens and particulate matter. Liquid or non-oil based particles from sprays that do not also emit oil aerosols or vapours

### Do not use for:

Do not use for gases and vapours, oil aerosols, asbestos, or sandblasting; particulate concentrations that exceed either 10 times the occupational exposure limit or applicable government regulations, whichever is lower.

### Warnings & Limitations

- 1) Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.
- 2) The respirator is designed to reduce exposures to certain airborne contaminants when fitted correctly. As such the respirator is not suitable for kids, those with beards or other facial hair or other conditions that prevent a good seal between the face and the sealing surface of the respirator.
- 3) The respirator does not eliminate the risk of disease or infection nor should it be used as a substitute for proper health management. Wearers should still take the necessary precautions when entering and working in a hazardous setting. Failure to follow all instructions and limitations can reduce the respirators effectiveness and may result in sickness or death
- 4) The respirator is designed as a single use disposable respirator - do not alter or wash this respirator.
- 5) Individuals with a compromised respiratory system, such as asthma or emphysema should consult a physician prior to use.

### Time Use Limitation


If respirator becomes damaged, soiled, or breathing becomes difficult, leave the contaminated area immediately and replace the respirator.

### User instructions

**Special or Critical User Instructions:** This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings as a Surgical Respirator conforming to recognized

standards for biocompatibility, flammability, and fluid resistance (at 120mmHg). This respirator has been tested to ASTM F-1862 to determine the resistance of this product to synthetic blood.

### Donning method

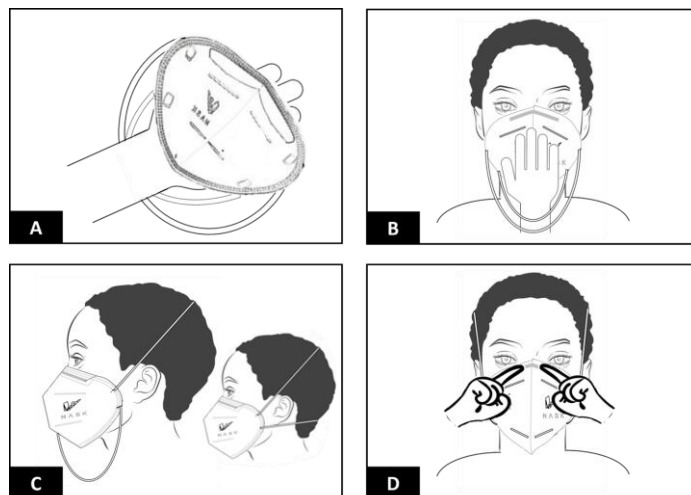

1. Hold the edges of the respirator with one hand, and position it over your nose and mouth
2. Use your other hand to pull the bottom headband over your head, positioning it around the base of your neck
3. Adjust the top headband over the crown of your head for a secure and comfortable fit.
4. With both hands, mold the metal nosepiece to the shape of your nose, ensuring a tight seal.
5. Perform a fit check by cupping your hands over the respirator and exhaling sharply. If you feel any air leakage around the edges, readjust the respirator and nosepiece for a better seal.
6. Avoid touching the respirator while wearing it.

### Doffing Method

1. Without touching the respirator and handling it only by the headband, gently lift the bottom headband from around your neck and pull it up and over your head.
2. Then lift the top headband and pull respirator
- Dispose of used respirator in accordance with applicable local regulation.
3. Perform hand hygiene