


FITTING INSTRUCTION

1. Hold the mask with nosepiece at fingertips.

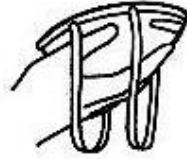

4. Place your fingertips from both hands at the top of the metal nosepiece.

2. Place the mask under your chin with nosepiece up.



5. Press the nosepiece to the shape of nose with two hands.

3. Pull the top strap over and rest at the top back of your head, the bottom strap around the neck below the ears.



6. To check the fit, place both hands over the mask and breathe vigorously. If air flows around the nose, readjust as in step 4 and fasten the straps.

Change respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted, or a proper face fit cannot be maintained. Remember, careful observance of these fitting instructions is an important step in safe respirator use.

To check the correct fit, cup with both hands the respirator and exhale powerfully. If any airflow indicates leakage around the edges of the respirator, stretch the headbands again by pulling them backwards. If necessary, the length of the elastic headband can be shortened by knotting. Then correct the position of the respirator (see picture) and check again.

WARNING! It is not allowed to use the mask without proper fit. If you can't achieve a proper seal, do not use the respirator.