USER'S INSTRUCTIONS

San-M Package Co., Ltd.  
N95 Respirator

TN02-29 (Medium)  
TN02-30 (Small)

San-M Package Co., Ltd.  
1086-1 Ojira, Shimada-city,  
Shizuoka 428-0009, Japan  
TEL : +81 547 45 4125

TN0229U1a  
[Rev. b] 10/16/2020

## USER'S INSTRUCTIONS

### WARNING!

Although this respirator is useful to protect the wearer from particulates, including those containing infectious disease, the risk for sickness or infection is not completely eliminated. Failure to read and follow all warnings and instruction will reduce respirator performance and may result in sickness and / or death. Before using, please read user's instruction.

### [USE INSTRUCTION]

This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.

Use for protection against particulate aerosols free of oil.

Do not use for gasses, vapors asbestos, paint spray or any oil-based aerosols, or in environments.

Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

Leave the use environment immediately in the case of dizziness, nausea, breathing difficulty or other symptoms of physical distress.

Don a new respirator if a proper fit cannot be obtained, if breathing becomes difficult, or the respirator becomes damaged.

Only for single use.

Store unused respirators in a clean, cool, dark, dry place.

## FITTING INSTRUCTION

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1. Separate mask to open.
2. Bend the nosepiece to conform mask to your nose.
3. Hold the center of headbands separately.
4. Hang the lower part on your chin and bring both bands over your head.
5. Place the upper band at your crown and the lower band at your neck.
6. Form the nosepiece tightly to make no chink.  
※ Using both hands, over the whole of the respirator to check the fit-breathing out to make sure that no breath escapes from around the nose, cheeks or edges of the respirator.

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Leave the use environment immediately in the case of dizziness, nausea, breathing difficulty or other symptoms of physical distress.


Don a new respirator if a proper fit cannot be obtained, if breathing becomes difficult, or the respirator becomes damaged.

Only for single use.

Store unused respirators in a clean, cool, dark, dry place.

## FITTING INSTRUCTION

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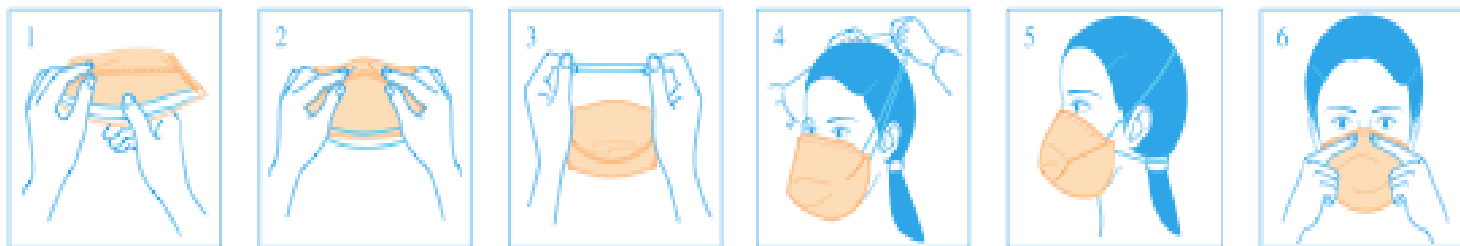1. Separate mask to open.
2. Bend the nosepiece to conform mask to your nose.
3. Hold the center of headbands separately.
4. Hang the lower part on your chin and bring both bands over your head.
5. Place the upper band at your crown and the lower band at your neck.
6. Form the nosepiece tightly to make no chink.

※ Using both hands, over the whole of the respirator to check the fit-breathing out to make sure that no breath escapes from around the nose, cheeks or edges of the respirator.

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