

RC-10

# Surgical N95 Particulate Respirator

TC-84A-9316



DCare AG  
Baarerstrasse 141, 6300 Zug,  
Switzerland  
E: info@rc-respirator.com  
T: +41 41 5014161



This respirator is approved only in the following configuration:

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-9316	Surgical/N95	RC-10 X	ABCJM NOS

### 1. Protection

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply

### 2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, PSHA and other applicable regulations
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in configuration as specified by the manufacturer.
- O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators
- S - Special or critical User Instructions and/or specific use limitations apply. Refer to User Instructions before donning.

S - Special or Critical User's Instructions: The model RC-10 filtering facepiece respirator has been manufactured by Approval Holder Makrite for DCare AG under TC-84A-9316.

S - Special or Critical User Instructions: This respirator has been approved as a NIOSH N95 filtering facepiece respirator, for use in healthcare settings, as a Surgical N95 Respirator conforming to recognized standards for biocompatibility, flammability, and fluid resistance at 160mmHg.

### SERVICE TIME LIMITATION:

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled or distorted. Do not use the same respirator again once it has been used in contaminated area.
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use, beginning from the first day of use)

### LIMITED WARRANTY


#### IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, expressed or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of DCare AG will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Dcare AG, products proved to be defectively manufactured. Except as provided above, DCare AG shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

RC-10      Made in China

## Instructions

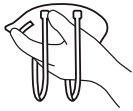

Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.



Press the respirator firmly against your face with the nosepiece on the bridge of your nose.



Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.



Using both hands, mould noseclip to the shape of the nose. Using one hand to adjust the nosepiece may result in improper fit.



To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit.

Caution: Respirator must be fit checked prior to each wearing.

### WARNINGS:

1. Wearer must read fitting instructions and use limitations before use.
2. This respirator is not guaranteed to protect against the risk of contracting a disease or infection.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

DCare AG  
Baarerstrasse 141, 6300 Zug,  
Switzerland  
E: info@rc-respirator.com  
T: +41 41 5014161