


JFY 44301 (P100)

Donning Procedures

1. Hold the respirator in your hand with the nosepiece at your fingertips and press it firmly against your face, with the nosepiece on the bridge of your nose.
2. With your other hand, take the lower strap, pull it over your head and position it below your ears.
3. Place the top strap high on the back of your head. Move the bottom strap over your head and position it below your ears.
4. Adjust the strap tension by pulling the straps as shown in the Figure.
5. Using both hands, mold the metal nosepiece to the shape of your nose and inhale vigorously. If no negative pressure is felt inside the respirator, adjust and tighten the nosepiece. If any leakage is detected, readjust the position of the respirator according to steps #4 and #5.

If you CANNOT achieve a proper fit, DO NOT enter the contaminated area. See your supervisor for proper fitting respiratory protection.

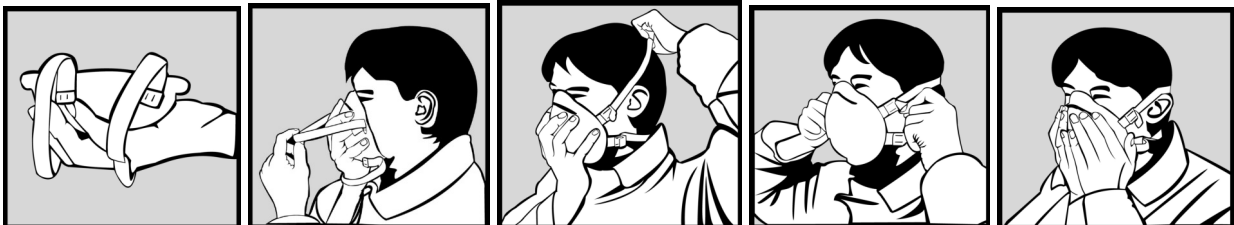

Gateway Safety 80904V (P100)

Donning Procedures

1. Hold mask in one hand with straps around back of hand. Place over mouth and nose with nose clip over bridge of nose. Pull lower strap over head and fit around neck. Pull upper strap over head and place it above ears.
2. Form nose clip around nose and adjust tension in straps by positioning straps on head until you obtain a positive face seal.
3. Perform a User Seal Check prior to each wearing. It is the wearer's responsibility to obtain a proper fit. Each time respirator is used, wearer must test fit by cupping both hands over the respirator and inhaling vigorously. If no negative pressure is felt inside the respirator, tighten the nosepiece and reposition or tighten straps until no leakage is detected. If you CANNOT achieve a proper face seal, DO NOT enter the contaminated area. See your supervisor.

If you CANNOT achieve a proper fit, DO NOT enter the contaminated area. See your supervisor for proper fitting respiratory protection.

