
INSTRUCTIONS / INSTRUCCIONES

To Put On This Respirator: Mise en place du respirateur: Para ponerse este respirador:

1. Hold the respirator in your hand with the metal adjusting strip away from you, allowing the headbands to hang below your hands.
2. Press the respirator firmly against your face with the narrow end on the bridge of your nose.
3. Pull the top strap over your head and position it high on the back of your head. Pull the bottom strap over your head and position it below your ears.
4. Using both hands, mold the metal strip to the shape of your nose.
5. Check the facepiece-to-face seal by cupping both hands over the respirator and inhale vigorously. A negative pressure should be detected inside the facepiece. If air leaks in from around your nose, reform the metal strip. If air leaks in from around the edges, reposition the straps. Repeat as necessary until no air leak is detected, and you have a satisfactory seal. If you cannot achieve a good face-to-facepiece seal do not enter the contaminated area and see your supervisor.

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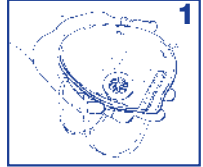Time Use Limitations:

The use and reuse of this respirator is subject to considerations of hygiene, damage and increased breathing resistance. Replace respirator if breathing becomes difficult. Service time should not be extended beyond 24 total hours of use.

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