



Shining Star Electronic Technology Co., Ltd.
9F, No. 33, Lane 166, Sec. 2, Meishi Rd., Yangmei Dist., Taoyuan, Taiwan 32668
TEL: +886 3 431 5486

This respirator is approved only in the following configuration:


TC-	Protection ¹	Respirator	Cautions and Limitations ²
		SS6001-N95	
84A-8125	N95	X	ABCJMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level)
Effective against particulate aerosols free of oil;
time use restrictions may apply

2. Cautions and Limitations

- A. Not for use in atmospheres containing less than 19.5% oxygen.
- B. Not for use in atmospheres immediately dangerous to life or health.
- C. Do Not exceed maximum use concentrations established by regulatory standards.
- J. Failure to properly use and maintain this product could result in injury or death.
- M. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O. Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators for use as surgical masks.

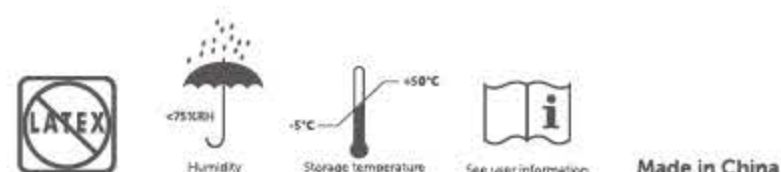

LIMITED WARRANTY IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face herof. The physical standards and specifications of Shining Star Electronic Technology will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Shining Star Electronic Technology, products proved to be defectively manufactured. Except as provided above, Shining Star Electronic Technology shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

⚠ WARNINGS:

1. Wearer must read fitting instructions and use Limitations before use.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

USER SEAL CHECK:



STEP 1.
Placed the respirator in hand with nosepiece facing downward, let the headband loose and hanging freely underneath your hand.



STEP 2.
Hold and press the respirator against your mouth, setting the nosepiece on the bridge of your nose.



STEP 3.
Pull and placed lower headband around your head and set it below your ears; pull and placed upper headband on the back of your head.



STEP 4.
Adjust the nosepiece and firmly press against your nose to form the shape to fit.



STEP 5.
To test fitness:
a. Press firmly with respirator and carefully adjust to center position against your mouth.
b. Inhale vigorously; if air is felt leaking into the respirator, readjust the headband or nosepiece for better fit. Try testing again after readjusting until a pass is achieved. Another method is to cover the facepiece with your hand and exhale vigorously, if air is felt leaking out through the edge of the respirator, readjust and try again.

Carefully follow the instruction for proper fitting to achieve maximum protection.

SERVICE TIME LIMITATION:

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, or distorted. Do not use the some respirator again once removed.
3. The maximum service time of this respirator is five (5) consecutive calendar day (including days of non-use), beginning from the first day of use.

Shining Star Electronic Technology Co.,Ltd.
9F, No. 33, Lane 166, Sec. 2, Meishi Rd., Yangmei Dist., Taoyuan, Taiwan 32668
TEL: +886 3 431 5486

Manufacturer : Hangzhou Ti Yun Industrial Co., Ltd.
First Building, Hang Xin Village, Canal Street, Yu Hang District, Hangzhou City, Zhejiang Province, 311103, China
TEL: +86 0571-8618-1168