


User's Instructions for Gateway Model No. 80202V

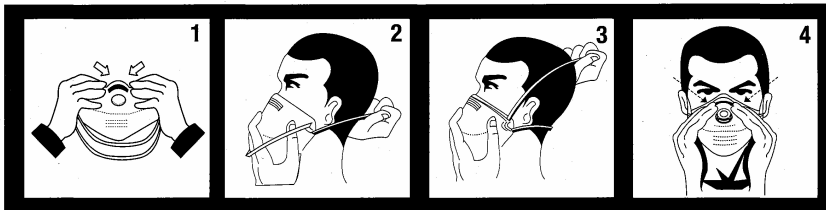

FITTING INSTRUCTIONS

1. Bend plastic nose clip to form nose shape.
2. Hold mask in one hand over mouth and nose with straps around back of hand. Pull lower strap over head and fit around neck. Pull upper strap over head and place it above ears.
3. Form nose clip around nose and adjust the tension in the upper and lower portions of the continuous strap to ensure a positive face seal.
4. Perform a User Seal Check prior to each wearing. It is the wearer's responsibility to obtain a proper fit. Each time respirator is used, wearer must test fit by cupping both hands over the respirator and exhaling vigorously. If air is felt flowing around nose or leaking around edges of the respirator, tighten the nosepiece and reposition or tighten strap until air is no longer escaping when wearer exhales. If you CANNOT achieve a proper face seal, DO NOT enter the contaminated area. See your supervisor.

WARNING

- A. Before occupational use of this respirator, a written respiratory protection program shall be implemented meeting all the requirements of OSHA 29 CFR 1910.134 such as training, fit testing, medical evaluation, and the applicable OSHA's substance specific standards. A properly selected respirator is essential to protect your health. Consult an Industrial Hygienist or Occupational Safety Professional to determine its suitability for your intended use.
- B. This respirator does not supply oxygen. Use only in well ventilated areas containing sufficient oxygen, at least 19.5%, to support life.
- C. Do not use when concentrations of contaminants are unknown, at levels dangerous to health and life, in excess of 10 times the permissible exposure limit (PEL) or in excess of other regulations.
- D. Leave contaminated area immediately if respirator is damaged, or if any physical distress such as breathing difficulty, dizziness or nausea occurs.
- E. Do not alter, modify or misuse respirator.
- F. Facial hair, beards or certain facial characteristics may prevent a proper fit and thus prevent the use of this respirator.
- G. Keep unused respirators in the display box, stored in a cool, dry, uncontaminated area, away from direct sunlight until use.
- H. Dispose of respirator after single use, increased breathing resistance or hygiene damage.

USE RESTRICTIONS

DO NOT USE in spray painting or sandblasting operations.

DO NOT USE for protection in asbestos environments.

DO NOT USE for protection against vapors, gases or oil-based aerosols.

DO NOT USE for welding applications or near any welding activities.

USE AGAINST non-oil particulate aerosols, such as coal, iron ore, cotton, flour, and for sprays that do not emit vapors or gases.

This Gateway Safety Model #80202V filtering face piece respirator is marketed and sold by Gateway Safety Inc. under a production agreement with San Huei United Co, Ltd. in accordance with NIOSH TC84A-3719.