



**aero  
pro™**

# AP0708B+ SURGICAL N95 RESPIRATOR FITTING INSTRUCTION

**1**




**Step 1:**  
Open the respirator.

**2**




**Step 2:**  
Pull both elastic straps  
out from the respirator.

**3**

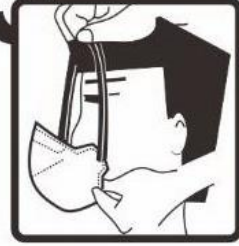

**Step 3:**  
Pass the elastics over  
the head and position  
the respirator under  
the chin with nose-clip  
facing upwards.

**4**

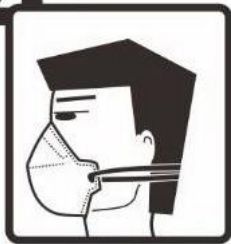

**Step 4:**  
Position the lower  
elastic strap around  
the neck and release it.

**5**



**Step 5:**  
Position the upper  
elastic strap on the top  
of the head and release  
it.

**6**



**Step 6:**  
Adjust the respirator  
over the face ensuring  
that pressing the nasal  
bar firmly and evenly  
against the nose for a  
snug and comfortable  
fit, so there are no  
leaks.

- Test the tightness of the respirator by exhaling or blowing out sharply. Check for leaks around the nose clip and the rest of the seal of the respirator. Any leaks around the seal must be eliminated before entering the contaminated area.
- It is the user's responsibility to ensure that the respirator provides the necessary level of protection for the type and concentration of the contaminants in the area where the respirator is intended for use.
- If you can't achieve a proper seal, do not use the respirator.