



AERO PRO CO., LTD.
No.90 , Hsin Kung 6 RD., Tien Chung Chen,
Chang Hua Hsien, Taiwan, R.O.C
886-4-875-6141



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	Protection ¹			Caution and Limitations ²
		AP112	AP112S	
84A-9726	N95	X	X	ABCJMNO


1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free oil; time use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5 percent oxygen.
- B - Not for use in the atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to User's Instruction, and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks

FITTING INSTRUCTION

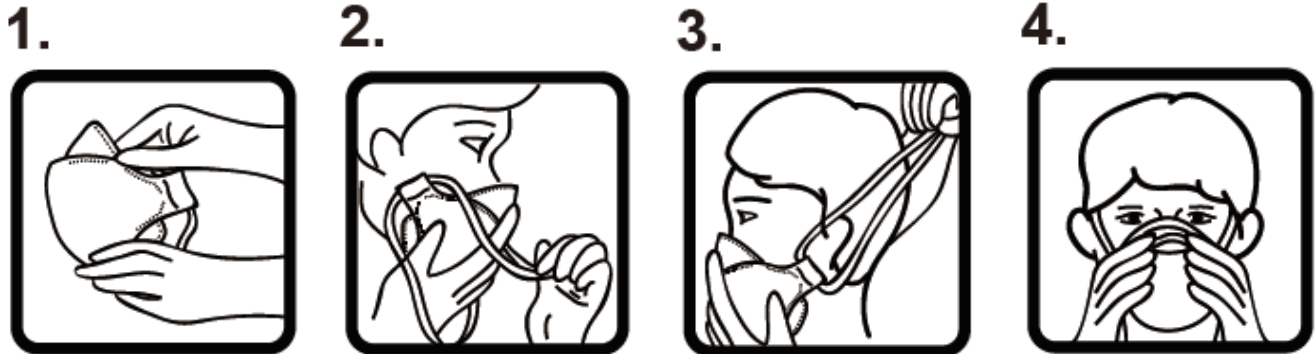

1. Make sure the respirator, including its head straps and nose clip, is not damaged or stained as it may negatively affect the efficiency of the respirator.
2. Place the respirator under your chin with the molded nose area up and pull the headbands over your head.
3. Position the lower head strap around the neck, below your ears, and the upper head strap around the back of your head.
4. Using both hands, adjust the nose clip to the shape of the bridge of the nose to ensure a close fit and good seal.

■ Change respirator immediately if breathing becomes difficult or respirator becomes damaged or distorted, or a proper face fit cannot be maintained. Remember, careful observance of these fitting instructions is an important step in safe respirator use.

■ To check the correct fit, cup with both hands the respirator and exhale powerfully. If any airflow indicates leakage around the edges of the respirator, stretch the headbands again by pulling them backwards. Then correct the position of the respirator and check again.

WARNING! If you can't achieve a proper seal, do not use the respirator.