PARTS LIST				REVISION
ITEM NO.	MATL. MASTER NO.	REQD.	DESCRIPTION	
				⑤ 03.21.06 A) ADDED INFLUENZA WARNING B) REVISED DRAWING LAYOUT DCN-1 ML BM

1. SIZE
  - 14-3/16" x 8-5/8"
  - DOUBLE FOLDED (BY VENDOR) TO 4-3/4" x 8-5/8 WITH COLORBANDS FACING OUT (FRONT) AND BARCODE FACING OUT (BACK)
  - SEE PAGE 4 FOR DIE LINE SPECIFICATIONS
2. MATERIAL
  - BOARD, WHITE, .014 SBS
  - UV COATING ONE SIDE (FRONT)
3. INK
  - ALL TEXT AND LINES PRINT PROCESS BLACK UNLESS OTHERWISE NOTED
  - ALL GRAPHICS PRINT STANDARD PANTONE PROCESS WAX FREE INK (CYAN, MAGENTA, YELLOW, BLACK) UNLESS OTHERWISE NOTED
  - GREEN PMS #3288 TO PRINT STANDARD PANTONE SPOT COLOR (NOT PROCESS PRINTING)
  - BLUE PMS #2945 TO PRINT STANDARD PANTONE SPOT COLOR (NOT PROCESS PRINTING)
  - INDUSTRIAL GRADE YELLOW PRINTS PMS #115
4. PRINTING
  - HEADER CARD FRONT PRINTS AS SHOWN ON PAGE 2
  - HEADER CARD BACK PRINTS INVERTED AS COMPARED TO PAGE 3
  - FINAL OUTCOME:
    - WHEN FOLDED HEADER CARD (WITH COLORBANDS FACING OUT) IS OPENED, INSIDE TEXT IS RIGHT READING
5. FOLDING
  - PANEL WITH NIOSH APPROVAL FOLDS FIRST; DOWN AND BEHIND BACK (BARCODE) PANEL
6. ARTWORK RETAINED ELECTRONICALLY BY MSA DESIGN RESOURCES
7. VENDORS AUTHORIZED IDENTIFICATION TO PRECEDE DATE REVISION CODE

AUTHORIZED COUNTRY	<b>MSA</b>	MINE SAFETY APPLIANCES CO. PITTSBURGH, PA 15230, U.S.A.	THIS DRAWING IS THE PROPERTY OF MINE SAFETY APPLIANCES COMPANY, PITTSBURGH, PENNSYLVANIA, AND IS MADE AVAILABLE TO YOU IN CONFIDENCE AND SUBJECT TO THE FOLLOWING: NO PERMISSION IS GRANTED TO PUBLISH, USE, REPRODUCE, TRANSMIT, OR DISCLOSE THIS DRAWING, OR ANY INFORMATION OR IDEA CONTAINED THEREIN TO OTHERS WITHOUT THE PRIOR WRITTEN CONSENT OF MINE SAFETY APPLIANCES COMPANY EXCEPT FOR THE MANUFACTURE OF ARTICLES FOR MINE SAFETY APPLIANCES COMPANY. SCHUTZVERMERK NACH DIN 34 BEACHTEN. MSAUER GmgH.		DO NOT SCALE DWG.
MSA STANDARD SHOP PRACTICES APPLY					
CODE SYMBOLS		JOB NO. E 002224	AGENCY APPROVAL CODE: NIOSH ATTR. DOC:		
CRITICAL	MINOR (NO MARK)	DRAWN BY ML 03.20.98	<b>CARD, HEADER, HARMFUL DUST DUST RESPIRATOR WITH EXHALATION VALVE, CONSUMER</b>		
MAJOR A		CHECKED BY TMS 04.03.98			
MAJOR B		PROJECT ENGINEER TMS 04.03.98			
100% INSPECTION		QUALITY ASSURANCE			
TOLERANCES UNLESS OTHERWISE SPECIFIED		MANUFACTURING	CADD REF: SCALE:		
FRACTIONAL ± 1/16					
2 PL. DEC. ± .06					
3 PL. DEC. ± .063					
ANGULAR ±					
SEE MSA STANDARD SHOP PRACTICES FOR DRILLED HOLE TOLERANCES		<b>C</b>	<b>817853</b>		SHEET 1 OF 4

Mine Safety Appliances Company  
Pittsburgh, Pennsylvania, USA  
1-800-672-2222

**AFFINITY PRO RESPIRATORS**

THESE RESPIRATORS ARE APPROVED  
ONLY IN THE FOLLOWING CONFIGURATIONS:

TC-	PROTECTION <sup>1</sup>	ALTERNATE AFFINITY PRO FILTERING RESPIRATOR		CAUTIONS AND LIMITATIONS <sup>2</sup>
		AFFINITY PRO N95 W/ VALVE STANDARD	AFFINITY PRO N95 W/ VALVE, SMALL	
		PR201-4	PR201-3	
84A-1018	N95	X	X	ABCJMNOP

**1. PROTECCION**

N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

**2. CAUTIONS AND LIMITATIONS**

- A - Not for use in atmospheres containing less than 19.5 percent oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configurations as specified by the manufacturer.
- O - Refer to User's Instructions, and / or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.

Mine Safety Appliances Company  
Pittsburgh, Pennsylvania, USA  
1-800-672-2222

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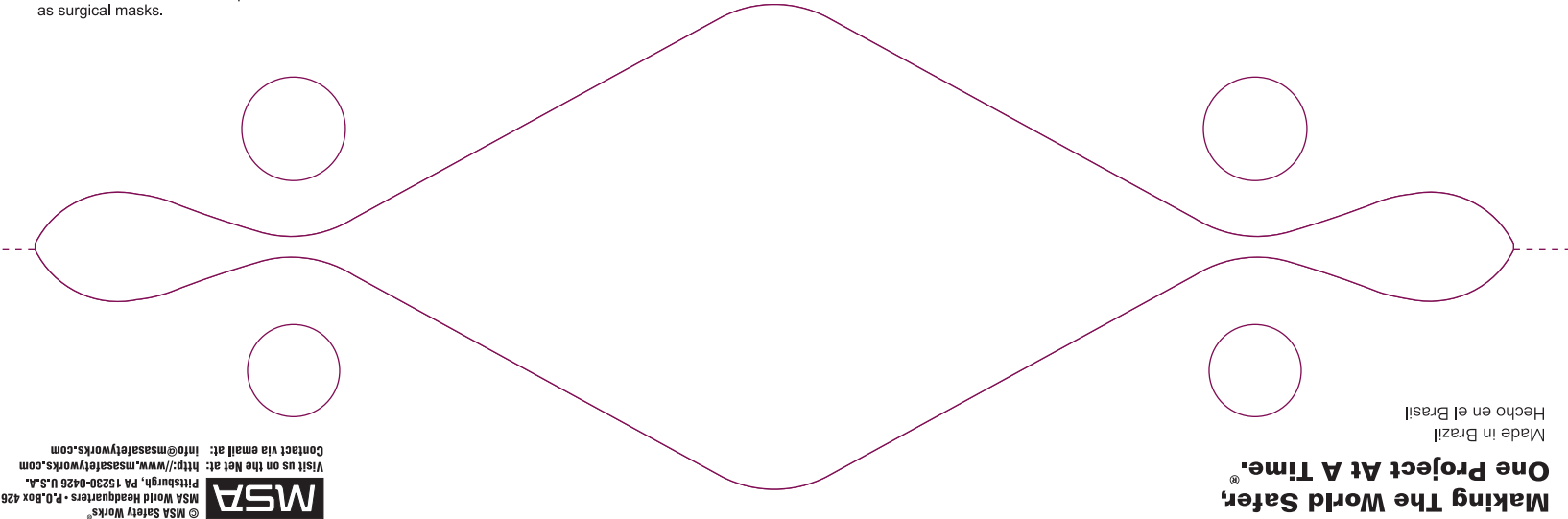
Making The World Safer.<sup>®</sup>  
One Project At A Time.  
Made in Brazil  
Hecho en el Brasil

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MSA warrants this product against defects in materials or workmanship during normal use for one year from date of purchase. Call 1 (888) 672-4692 and MSA will repair or replace the product, at MSA's sole discretion. MSA shall not be liable for incidental or consequential damages resulting from the use of this product. MSA specifically disclaims any warranty of merchantability or fitness for a particular purpose. Some states do not allow limitations on how long an implied warranty lasts or the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may also have other rights which vary from state to state.

**WARNING / ADVERTENCIA**  
DO NOT use while performing or if near sandblasting, concrete or masonry operations, fire fighting, or any operation that may produce airborne particles containing silica (sand), asbestos, lead, or other substances reported to cause cancer. Read all instructions and warnings before using. Misuse can result in serious personal injury, illness or death. See inside for NIOSH approvals.  
NO utilice este producto al realizar o al estar cerca de procedimientos donde se utilice chorro de arena, hormigón o albañilería, al combatir incendios o durante cualquier otro procedimiento que genere lesiones personales graves, enfermedades o incluso la muerte. Vea adentro para las aprobaciones de NIOSH.

**Questions? Comments? Call toll-free 1-888-MSA-HOW2 (1-888-672-4692).**  
**Si tiene preguntas o comentarios, llame si 1-888-MSA-HOW2 (1-888-672-4692).**

**SAFETY WORKS**  
SINCE 1914  
Dust Respirator with exhalation valve  
Model 00817626


# Dust Respirator

with exhalation valve

RECOMMENDED BY MSA **for Harmful Dust**

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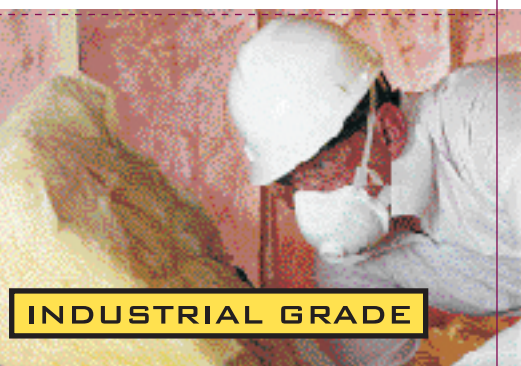**INDUSTRIAL GRADE**

RESPIRATORY PROTECTION  
**Sanding/Saw Dust, Fiberglass, Drywall Dust**

**Uses:**

- Sanding/sawing wood
- Installing fiberglass insulation
- Cutting or sanding drywall

**Meets OSHA and NIOSH Requirements**

**Benefits:**

- Exhalation valve allows cooler, easier breathing
- 2 adjustable straps and patented shape provide better seal
- Foam sealing surface improves comfort/fit
- Reclosable package for storage

Contents/Contenido: 1

00817626

## INSTRUCTIONS

### IMPORTANT

Fitting instructions to be followed each time the respirator is worn:

1. Hold the respirator in your hand, allowing the headbands to hang freely below your hands.
2. Position the respirator comfortably under your chin and against the bridge of your nose.
3. Pull the top headband over your head, placing it at the crown of your head. Pull the bottom headband over your head and place it below your ears.
4. FIT CHECK: Place both hands completely over the respirator. Inhale sharply. Be careful not to disturb the position of the respirator. A negative pressure (mask wanting to collapse) should be felt. If not, or if air leaks around the edges of the respirator, reposition the respirator and headbands for a better fit. If unable to adequately seal filter with hands, obtain the optional respirator fit check cup to perform this test.

If you CANNOT achieve a proper fit, DO NOT enter the contaminated area. See your supervisor or call 1-888-MSA-HOW2 (1-888-672-4692) for proper fitting respiratory protection.

### Applications

Particulate aerosols free of oils meeting the requirements described herein.

### Preparations for Use:

1. A properly selected respirator is essential to protect your health. Before using this product, consult your supervisor and review the instructions. If uncertain about the suitability of this respirator for your intended use, call MSA Safety Works at 1-888-MSA-HOW2 (1-888-672-4692).
2. Before use, the wearer must be trained by the employer in proper respirator use.
3. OSHA (29 CFR 1910.134) requires that each wearer of this respirator have a qualitative or quantitative fit test performed to determine the amount of protection it will provide.

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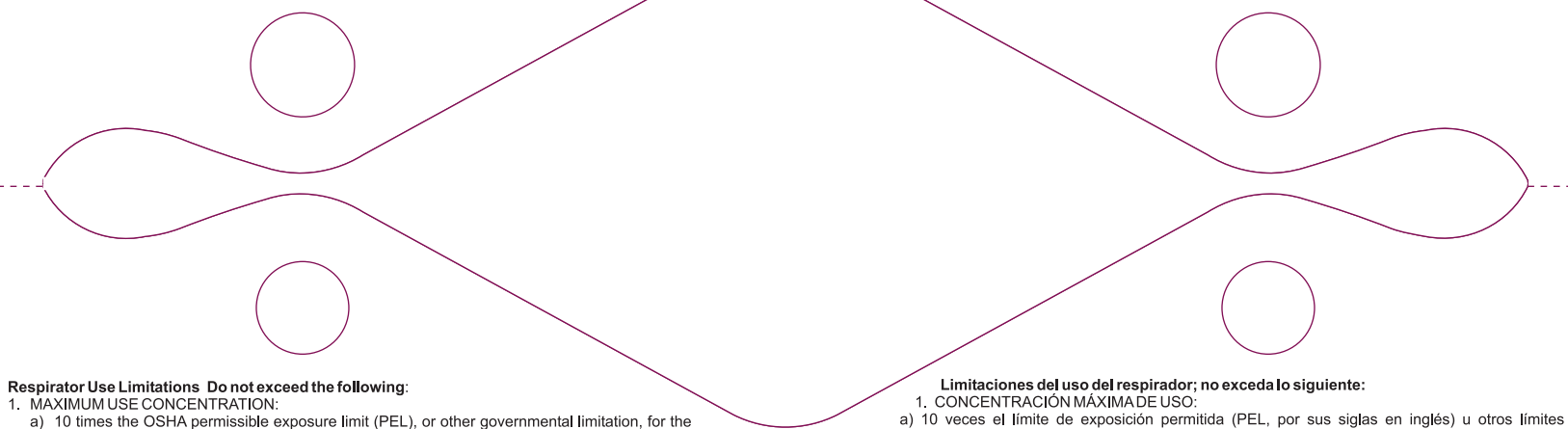
### Respirator Use Limitations Do not exceed the following:

1. MAXIMUM USE CONCENTRATION:
  - a) 10 times the OSHA permissible exposure limit (PEL), or other governmental limitation, for the contaminant(s) present.
  - b) Immediately dangerous to life or health (IDLH) concentration for any contaminant present.
2. 8 hours of use (continuous or intermittent). [Service time can be extended by performing an evaluation in the specific workplace setting that demonstrates (a) that the extended use will not degrade the filter efficiency below 95%, or (b) that the total mass loading of the filter is less than 200 mg.]
3. Any applicable regulatory or industry standard limitation or any limitation contained in the above warning.

### Replace the respirator when:

- Unable to pass a fit check.
- Breathing resistance is too high.
- Time use limitation is reached.

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### ⚠ WARNING

This respirator provides LIMITED protection. It may help reduce exposure to airborne biological agents, including avian (bird) flu virus, other types of influenza, SARS, or other bacterial or viral biological agents, but will not eliminate the risk of exposure, infection, illness, or death. • This respirator can filter these contaminants at the designated efficiency; however, the government has NOT established a safe level of exposure to these substances. Therefore, breathing a particle of any of these substances may be unsafe. • The Center for Disease Control and Prevention (CDC) recommends health-care workers exposed to persons under evaluation for avian influenza use a respirator at least as protective as a NIOSH approved N-95 respirator. They specify that respirators MUST be used in accordance with a complete respiratory protection program as required by OSHA, which includes proper selection, training, fit-testing, and fit-checking. Detailed information on a respiratory protection program is available by contacting OSHA or visiting [www.osha.gov](http://www.osha.gov). • Do NOT remove respirator in contaminated areas. The outer surface of the respirator MUST be treated as if it is contaminated at all times. A respirator will NOT prevent one from catching the flu in other ways such as by touching the mouth, nose, or eyes with contaminated hands or objects. Biological agents, such as the flu virus, can be transmitted when infected individuals cough or sneeze and spread virus particles through the air to exposed surfaces which are touched. Tight-fitting safety goggles, or a full-facepiece respirator, MUST be used to help prevent transmission of the disease. • The CDC recommends frequent hand washing and wearing gloves to help prevent transmission of the disease due to exposure to surfaces where contaminants may be present, and also immediately following removal of the respirator. • Do NOT reuse or share maintenance-free respirators. ALWAYS clean cartridge-style respirators before reuse in accordance with the instructions provided. • This respirator was NOT evaluated for use by children or by people with a medical condition that might prevent the use of a respirator. These individuals MUST consult a physician before using. • Failure to follow all warnings and instructions can result in serious personal injury or death.

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### ⚠ WARNING

1. Failure to select the appropriate respirator for all contaminants and their concentrations, or failure to follow all warnings and instructions including the accompanying NIOSH Cautions and Limitations, may result in exposure to hazardous materials and risk serious injury, illness or death.
2. DO NOT use under any of the following conditions:
  - a) For protection against silica (sand), asbestos, lead, or other substances reported to cause cancer.
  - b) For protection against gases or vapors, paint spray, unknown contaminants, or any contaminant whose concentration is unknown or immediately dangerous to life or health (IDLH).
  - c) In areas not adequately ventilated or containing less than 19.5 percent oxygen.
  - d) While performing or if near sandblasting, concrete or masonry operations, building demolition, firefighting, or any operation that may produce airborne particles containing substances reported to cause cancer. Small particles are not visible and may remain airborne even after operation has ceased. DO NOT use where these airborne particles may be present.
  - e) In atmospheres where the concentrations of contaminants exceed the Maximum Use Concentration or if the governmental exposure limit (PEL, REL, TLV) is not known.
  - f) When conditions prevent a good face seal such as beards, large sideburns, deep scars, missing teeth, or other facial characteristics that prevent direct contact between the skin and the sealing surface of the respirator.
  - g) If not trained in respirator use.
3. Immediately leave contaminated area if:
  - Breathing becomes difficult;
  - Dizziness or other distress occurs;
  - You smell, taste, or sense irritation from the contaminant;
  - The respirator becomes damaged;
  - You develop an allergic reaction to the respirator material.
4. This respirator will NOT eliminate exposure to contaminants. It is designed to reduce the breathing of contaminants by the wearer to an acceptable concentration when properly selected, fitted, and used. If a higher level of respiratory protection is desired, use an air-supplied respirator or full facepiece air-purifying respirator.
5. NEVER alter or modify this respirator.

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