


NIOSH N95 FLAT FOLDED RESPIRATOR

By **AllergyZone**

FITTING INSTRUCTIONS:

1

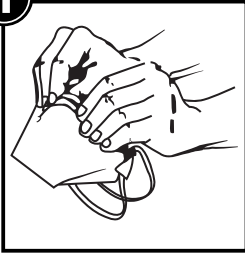

Step 1

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

2



Step 2

Open the mask. Let your chin into the bottom of the mask. Stretch and position the lower strap over your ears and upper strap over the head. Using both hands, mold the metal nosepiece to the shape of your nose.

3



Step 3

To test fit: (a) cup both hands over the mask being careful not to disturb position, (b) inhale vigorously, (c) exhale vigorously. If air leaks around the edges during inhalation or exhalation, re-position the straps or adjust strap tension to better fit. Repeat (a), (b) and (c). If there are still noticeable leaks then this mask may not be suitable for you and we advise you to seek another style from our range.

4



Step 4

Careful observance of these fitting instructions is an important step in safe respirator use.