


Fitting Instruction

* Donning procedure



1. Please inspect the respirator prior to each use to assure there is no damage, and then separate respirator to open.

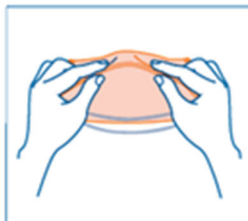

2. Bend the nosepiece to conform respirator to your nose.



3. Hold the center of headbands separately.



4. Hang the lower part on your chin and bring both bands over your head.



5. Place the upper band at your crown and the lower band at your neck.



6. Form the nosepiece tightly to your face by using the fingers of both hands, form the nose piece to the nose.

※ To check respirator fit: Cup both hands around the respirator, being careful not to disturb the position. Breathe out, if no breath escapes from around the nose, cheeks or edges of the respirator the respirator is fitted properly. If air leaks around the facepiece, reposition the straps and try again. Follow these fitting instructions during each use to achieve a proper fit.

* Doffing procedure



1. Grab both elastic band with both hands, and pull it backwards to remove it from the head.