

140 MM

Respiratory Protection



130 MM

140 MM

130 MM

Jet Black
Pantone
Process Black C
CMYK 0 0 0 100
RGB 0 0 0
HEX #000000

Ferrari
Pantone
185 C
CMYK 0 100 80 0
RGB 237 25 58
HEX #ED193A

Pearl
Pantone
White
CMYK 0 0 0 0
RGB 255 255 255
HEX #ffffff

WORKSafe®
PEOPLE PROTECTION

WORKSafe®
PEOPLE PROTECTION

N95 Respirators



WSRM211N95

contain 20
respirators



WSRM111N95
(small)

contain 20
respirators



WSRM512N95

contain 10
respirators



WSRM248N95

contain 10
respirators



PDS International Pte Ltd
10 Pandan Crescent #05-03/04 (LL2)
Singapore 128466
+65 6776 6200



THIS RESPIRATOR IS APPROVED IN THE FOLLOWING CONFIGURATIONS:

TC-	Protection ¹	Respirator				Cautions & Limitations ²
		WSRM211N95	WSRM111N95	WSRM512N95	WSRM248N95	
84A-7867	N95	X				ABCJMNOPS
84A-7871	N95		X			ABCJMNOPS
84A-7869	N95			X		ABCJMNOPS
84A-7882	N95				X	ABCJMNOPS

- 1 PROTECTION**
N95 - Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; Time use restrictions may apply
- 2 PLEASE TAKE NOTE**
- A Not for use in atmospheres containing less than 19.5 percent oxygen
 - B Not for use in atmosphere immediately dangerous to life or health
 - C Do not exceed maximum use concentrations established by regulatory standards
 - J Failure to properly use and maintain this product could result in injury or death
 - M All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulation
 - N Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer
 - O Refer to user instructions, and/or maintenance manuals for information on use and maintenance of these respirators
 - P NIOSH does not evaluate respirators for use as surgical masks
 - S Special or critical User's Instructions and/or specific use limitation apply. Refer to User's Instructions before donning

WARNINGS:
1) This respirator should not be used to protect against chemical warfare agents or nuclear material
2) Please read fitting instructions and usage limitations before using this product
3) Before use, the wearer must be trained in the proper use, including fit testing, in accordance with applicable safety and health standards for the containment and exposure level in the work area as per the OSHA guidelines, 29 CFR 1910.134.

worksafe.com.sg

WORKSafe®
PEOPLE PROTECTION

WORKSafe®
PEOPLE PROTECTION

N95 Respirators



WSRM211N95

contain 20
respirators



WSRM111N95
(small)

contain 20
respirators



WSRM512N95


contain 10
respirators



WSRM248N95

contain 10
respirators

FITTING INSTRUCTIONS

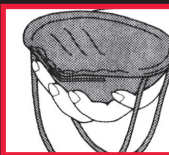

1 Cup the respirator face-down in your hand.



2 Place the respirator under your chin. The top strap goes around the upper back portion of your head, while the lower strap goes under your neck and below your ears.



3 Use both hands to mold the nose area into the shape of your nose. Push inwards while contouring the sides of the nosepiece.



4 Check the seal by placing both hands over the respirator and exhale. If air leaks around the nose, redo Step 3. If air leaks around the edges, ensure the straps are secure.



Service Time Limitation:

- When in dirty conditions, do not use the same respirator for more than 8 hours in total.
- If the respirator or your breathing is compromised, leave the area and remove the respirator.
- The maximum service time is 5 consecutive calendar days, beginning on the first day of use.

Limited Warranty:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of PDS INTERNATIONAL PTE LTD will be met by products.

Exclusive Remedies:

Damages for the breach of this limited warranty are limited to the replacement of such quantity of PDS INTERNATIONAL PTE LTD products proven to be defectively manufactured. Except as provided above, PDS INTERNATIONAL PTE LTD shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

S - Special or Critical User's Instructions: The model WSRM211N95 & WSRM111N95 & WSRM512N95 & WSRM248N95 filtering facepiece respirator has been manufactured by Approval Holder Makrite for Company PDS INTERNATIONAL PTE LTD under TC-84A-7867 & TC-84A-7871 & TC-84A-7869 & TC-84A-7882.

If you **CANNOT** achieve proper seal, DO **NOT** enter the isolation or treatment area. See your supervisor.

WORKSafe®
PEOPLE PROTECTION

LOT NO.