

USER'S INSTRUCTIONS

Unicharm Corporation
N95 Respirator

SNF-M01 Medium
SNF-S01 Small

TN0227PL2UIb

[Rev.b] 10/06/2023

Unicharm Corporation
3-5-19, Mita, Minato-ku,
Tokyo, 108-8575 Japan

TEL: + 81-3-3451-5111

These SNF-M01 Medium and SNF-S01 Small filtering face piece respirators” have been manufactured by Approval Holder San-M Package Co., Ltd. for Unicharm Corporation under TC-84A-9252.

USER'S INSTRUCTIONS

WARNING!

Although this respirator is useful to protect the wearer from particulates, including those containing infectious disease, the risk for sickness or infection is not completely eliminated. Failure to read and follow all warnings and instruction will reduce respirator performance and may result in sickness and / or death. Before using, please read user's instruction.

[USE INSTRUCTION]

This respirator does not supply oxygen and must not be used in atmospheres containing less than 19.5% oxygen.

Use for protection against particulate aerosols free of oil.

Do not use for gasses, vapors asbestos, paint spray or any oil-based aerosols, or in environments.

Before occupational use of this respirator a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

Leave the use environment immediately in the case of dizziness, nausea, breathing difficulty or other symptoms of physical distress.

Don a new respirator if a proper fit cannot be obtained, if breathing becomes difficult, or the respirator becomes damaged.


Only for single use.

Store unused respirators in a clean, cool, dark, dry place.


S-[SPECIAL OR CRITICAL USER'S INSTRUCTIONS]

Fitting Instruction


* Donning procedure

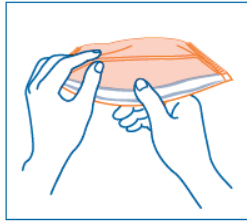

1. Please inspect the respirator prior to each use to assure there is no damage, and then separate respirator to open.

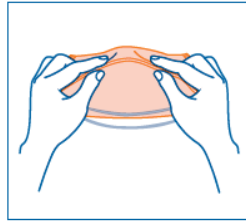


2. Bend the nosepiece to conform respirator to your nose.

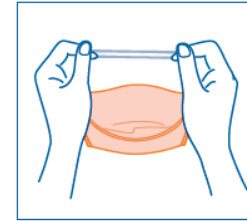

3. Hold the center of headbands separately.



4. Hang the lower part on your chin and bring both bands over your head.

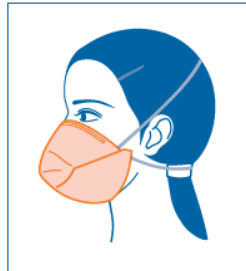

5. Place the upper band at your crown and the lower band at your neck.



6. Form the nosepiece tightly to your face by using the fingers of both hands, form the nosepiece to the nose.

※ To check respirator fit: Cup both hands around the respirator, being careful not to disturb the position. Breathe out, if no breath escapes from around the nose, cheeks or edges of the respirator the respirator is fitted properly. If air leaks around the facepiece, reposition the straps and try again. Follow these fitting instructions during each use to achieve a proper fit.

* Doffing procedure



1. Grab both elastic band with both hands, and pull it backwards to remove it from the head.