


# 8


## STEPS

TO  
DON  
THE

## SHIELD-95-PRO RESPIRATOR

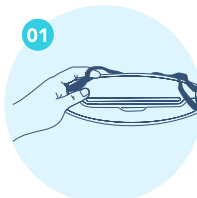

Turn the respirator so that you are looking at the back side.

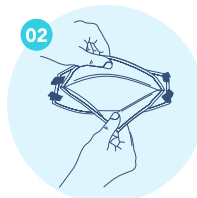

Pull up on the nose piece and down on the chin pieces to create an opening.



Place the respirator on your face so that the center of the nose piece wire is against the bridge of your nose and the chin piece is below your chin.



Pull both straps over your head. Move the lower strap around your neck; pull the top strap onto your head.



Using one hand, pinch the nose wire across the bridge of your nose and along the sides.



Using both hands, press each side of the nose wire strip to conform to your cheek bones on both sides.



With the respirator fitted to your face, place finger and thumb on the ends of the ClearView Bar and pinch towards the face to increase downward visibility.



Seal Check: Perform a seal check by covering the front of the respirator with both hands, being careful not to disturb the position of the respirator. Exhale sharply.

