

尺寸:140X130X200MM
颜色: CMYK

140 MM



130 MM

140 MM

130 MM

200 MM

Particulate Respirator

N95

Better Facemask... Powerful Protection!




Great Hawk™

20/Box

Read Instructions Before Use

MANUFACTURE DATE
EXPIRY DATE
LOT#
NIOSH LOT#



Great Hawk Air Service

A subsidiary of American Kono Group

701 5th Ave, FL 42

Seattle, WA 98104


USA

Fax: 1-360-368-5999

Email: info@N95-US.com

Website: www.N95-US.com

Page: www.facebook.com/GHA800



THIS RESPIRATOR IS APPROVED ONLY IN THE FOLLOWING CONFIGURATION:

TC-	PROTECTION¹	RESPIRATOR	CAUTIONS AND LIMITATIONS²
		GH9500	
84A-5411	N95	X	ABCJMNOPS

1.Protection

The NIOSH-approved N95 Respirator effectively provides the wearer with at least 95% filter efficiency level against airborne particles as small as 0.3 micron free of oil. Time use restrictions may apply.

2.Cautions and Limitations

A. Not for use in atmospheres containing less than 19.5 percent oxygen.
B. Not for use in atmospheres immediately dangerous to life or health.
C. Do not exceed maximum use concentrations established by regulatory standards.
J. Failure to properly use and maintain this product could result in injury or death.
M. All approved respirators shall be selected, fitted, used, and maintained in accordance with NIOSH, OSHA and other applicable regulations.
N. Never substitute, modify, and/or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O. Refer to user's instructions and/or maintenance manuals for information on use and maintenance of these respirators.
P. NIOSH does not evaluate respirators for use as surgical masks.
S. Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S. Special or Critical User's Instructions: The model GH9500 N95 filtering facepiece has been manufactured by Makrite for private label holder Great Hawk Air Service under TC-84A-5411. Great Hawk Air Service is a subsidiary of American Kono Group, USA.

SERVICE TIME LIMITATION

- Do not use the same respirator for more than 8 hours, continuous or intermittent in a single day, in dirty workplaces that could result in high filter loading.
- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once it has been used in contaminated area.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

Particulate Respirator

N95

Better Facemask... Powerful Protection!



Great Hawk™

IMPORTANT NOTICE TO BUYER'S AND END-USER'S

The Manufacturer warrants this personal protective equipment ("PPE") to be free from defect/s in material and workmanship.

In the event that this PPE is proven, upon Manufacturer's inspection, to contain defect/s in either material or workmanship, the Manufacturer may, at its option, to either replace any portion proving defective or to refund to the direct Buyer the original purchase price paid by the direct Buyer for this PPE.

This warranty and Buyer's remedy of replacement or refund constitute the Manufacturer's exclusive warranty and the Buyer's sole remedy, in lieu of any other warranties expressed or implied.

Except as provided above, the Manufacturer shall not be liable for any loss or damage, whether direct, indirect, consequential, or incidental occasioned by the sale, use or misuse of this PPE, other than death or personal injury directly resulting from the Manufacturer's negligence.


Buyer shall determine the suitability of this PPE for its intended use, and except provided above, assumes all risks and liability in connection herewith.

WARNINGS:
1. Wearer must read fitting instructions and use limitations before use.
2. This respirator is not guaranteed to protect against the risk of contracting a disease or infection.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing.

Made in China



Fitting Instructions

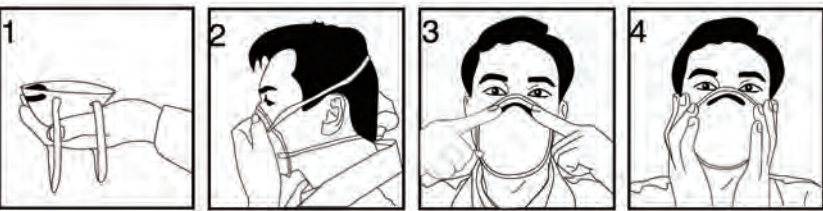

- Hold respirator in hand with nose clip at your fingertips and let the straps hang downwards below hand.
- Place the respirator under the chin and its nose clip upwards. Pull the upper strap to the top back of your head and then pull the shorter lower strap to the back of neck, below the ears.
- Push the nosepiece inwards with your fingertips to conform the shape of your nose.
- Face fit check must be done each time the respirator is worn as follows:
 - Cover the respirator with both hands, breathing sharply. If leakage along the nose bridge is detected, repeat step 3 and adjust the shape of clip to the nose bridge.
 - If leakage is detected in the rest of respirator edge, adjust the straps again.
 - If leakage cannot be eliminated, do not wear the respirator and stop going into the contaminated environment and turn to the supervisor for help.

Great Hawk Air Service
A subsidiary of American Kono Group
701 5th Ave, FL 42
Seattle, WA 98104
USA
Fax: 1-360-368-5999
Email: info@N95-US.com
Website: www.N95-US.com
Page: www.facebook.com/GHA800