
FITTING INSTRUCTIONS

Must be followed each time respirator is worn.

1. Hold the respirator in your hand with the metal adjusting strip away from you, allowing the headbands to hang below your hands.
2. Press the respirator firmly against your face with the narrow end on the bridge of your nose.
3. Pull the top strap over your head and position it high on the back of your head. Pull the bottom strap over your head and position it below your ears.
4. Using both hands, mold the metal strip to the shape of your nose.
5. Check the facepiece-to-face seal by cupping both hands over the respirator and inhale vigorously. A negative pressure should be detected inside the facepiece. If air leaks in from around your nose, reform the metal strip. If air leaks in from around the edges, reposition the straps. Repeat as necessary until no air leak is detected, and you have a satisfactory seal. If you cannot achieve a good face-to-facepiece seal do not enter the contaminated area and see your supervisor.

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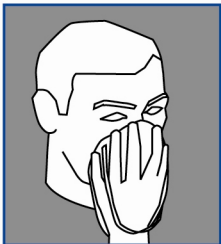Time Use Limitations:

The use and reuse of this respirator is subject to considerations of hygiene, damage and increased breathing resistance. For dirty workplaces, service time should not be extended beyond 8 total hours of use.

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Special Instructions:

The 1111 N95 Respirator has been manufactured by Aswan for OKI Bering under TC-84A-4659. NOTE: Before issuing any respirator to be worn, a face fit test must be performed either qualitatively or quantitatively. See OSHA regulations 1910.134. Mask should not be worn for more than 8 hours.

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