

140 MM

BHID-S302

20

Blancos • Con Elástico

White • with Elastic

N95 Particulate Respirators  
N95 Respiradores para partículas

BLACKHAWK  
INDUSTRIAL  
BIG ENOUGH TO SERVE | SMALL ENOUGH TO CARE

BLACKHAWK  
INDUSTRIAL  
BIG ENOUGH TO SERVE | SMALL ENOUGH TO CARE



N95 Particulate Respirators  
N95 Respirador para partículas

White • with Elastic

Blancos • Con Elástico



739630

901363

REORDER# /  
NÚMERO DE ORDEN  
BHID-S302

20

MASKS/MÁSCARAS

Lot Number/MFG Date:  
Número de Lote/Fecha de Fabricación:

Expiration Date:  
Fecha de Caducidad:

130 MM

140 MM

130 MM

200 MM



BlackHawk Industrial  
Broken Arrow, OK 74012  
1-855-610-1001  
www.BHID.com

NIOSH

This respirator is approved only in the following configuration:

TC*	Protection <sup>1</sup>	Respirator	Cautions & Limitations <sup>2</sup>
84A-5411	N95	BHID-S302 X	ABCJMNOPS

**1 Protection:**

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

**2 Cautions & Limitations:**


- A - Not for use in atmospheres containing less than 19.5% oxygen.
- B - Not for use in atmospheres immediately dangerous to life or health.
- C - Do Not exceed maximum use concentrations established by regulatory standards.
- J - Failure to properly use and maintain this product could result in injury or death.
- M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
- N - Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O - Refer to users' instructions and/or maintenance manuals for information on use and maintenance of these respirators.
- P - NIOSH does not evaluate respirators for use as surgical masks.
- S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model BHID-S302 filtering face piece respirator has been manufactured by Approval Holder Makrite for BlackHawk Industrial under TC-84A-5411.

**SERVICE TIME LIMITATION:**

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained, or if the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

Distributed by: BlackHawk Industrial  
Broken Arrow, OK 74012  
1-855-610-1001 • www.BHID.com

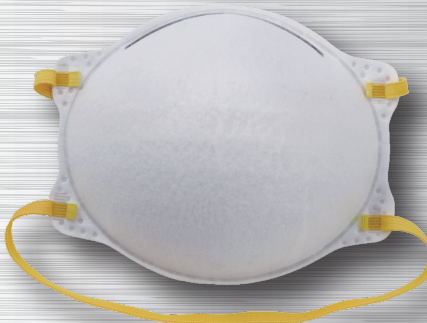

N95 Particulate Respirators  
N95 Respirador para partículas

White • with Elastic

Blancos • Con Elástico

**LIMITED WARRANTY | IMPORTANT NOTICE TO THE PURCHASER:**

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of BlackHawk Industrial will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of BlackHawk Industrial products proved to be defectively manufactured. Except as provided above, BlackHawk Industrial shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.

**WARNINGS:**

1. Wearer must read fitting instructions and use limitations before use.
2. This respirator does not protect against the risk of contracting disease or infection.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134, which includes medical evaluation, training and fit testing.

Made in China  
Hecho En China



**RESPIRATOR FITTING INSTRUCTIONS:**

1



**STEP 1:**

Hold the respirator in hand with the nose piece at your fingertips, allowing the headbands to hang freely below your hand.

2



**STEP 2:**

Press the respirator firmly against your face with the nose piece on the bridge of your nose.

3



**STEP 3:**

Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.

4



**STEP 4:**

Using both hands, mold nose piece to the shape of your nose.

5



**STEP 5:**

To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit.

Please carefully follow these fitting instructions during each use to achieve proper fit.