




SIZE:140*130*175mm



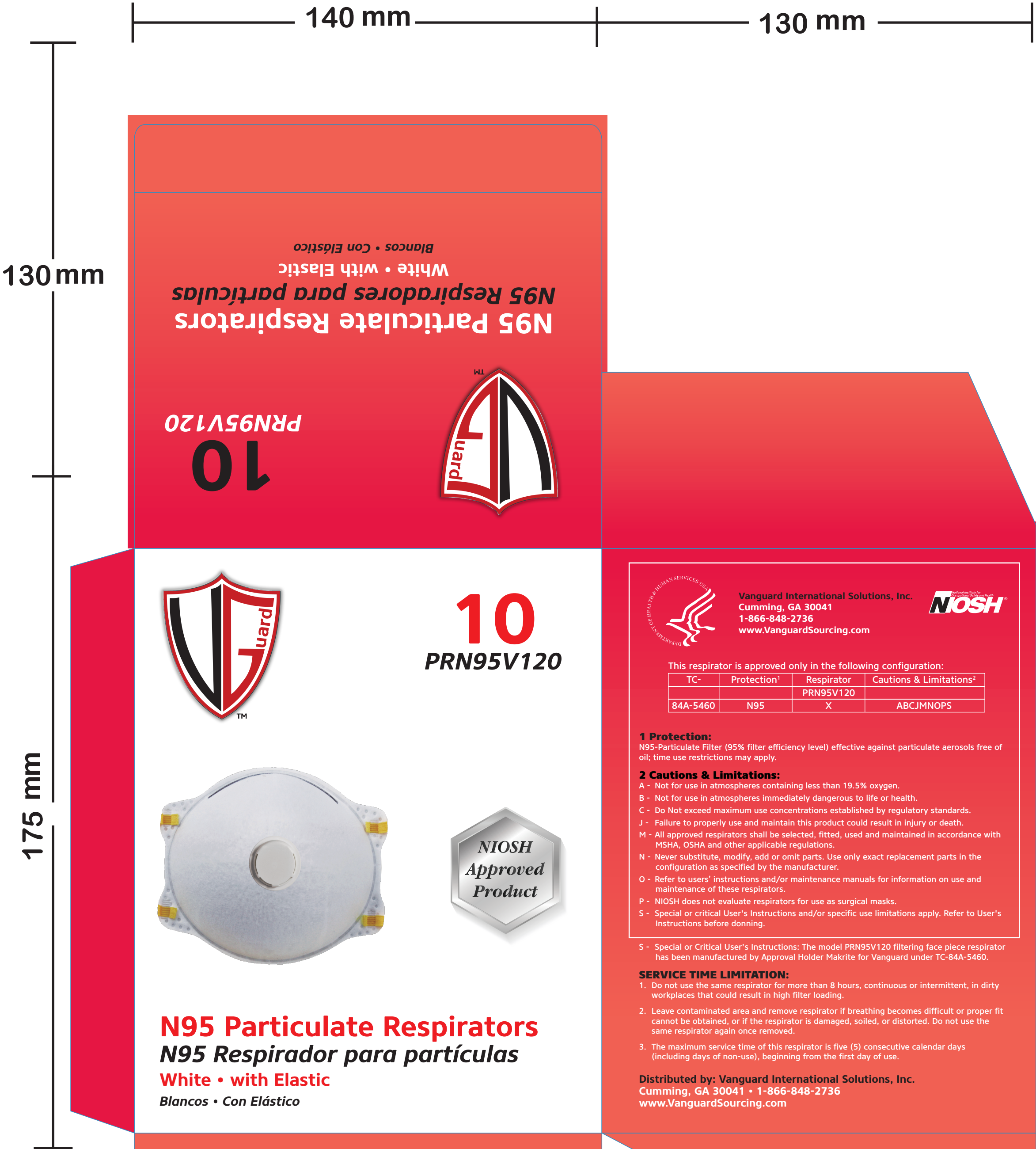

Vanguard International Solutions, Inc.
Cumming, GA 30041
1-866-848-2736
www.VanguardSourcing.com



This respirator is approved only in the following configuration:

TC-	Protection¹	Respirator	Cautions & Limitations²
		PRN95V120	
84A-5460	N95	X	ABCJMNOPS

1 Protection:
N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

2 Cautions & Limitations:
A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do Not exceed maximum use concentrations established by regulatory standards.
J - Failure to properly use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
N - Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to users' instructions and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.
S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model PRN95V120 filtering face piece respirator has been manufactured by Approval Holder Makrite for Vanguard under TC-84A-5460.

SERVICE TIME LIMITATION:

- Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained, or if the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

Distributed by: Vanguard International Solutions, Inc.
Cumming, GA 30041 • 1-866-848-2736
www.VanguardSourcing.com



10
PRN95V120

N95 Particulate Respirators
N95 Respirador para partículas
White • with Elastic
Blancos • Con Elástico

LIMITED WARRANTY | IMPORTANT NOTICE TO THE PURCHASER:
This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Vanguard International Solutions, Inc. will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Vanguard International Solutions, Inc. products proved to be defectively manufactured. Except as provided above, Vanguard International Solutions, Inc. shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.

WARNINGS:

- Wearer must read fitting instructions and use limitations before use.
- This respirator does not protect against the risk of contracting disease or infection.
- Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134, which includes medical evaluation, training and fit testing.

Made in China
Hecho En China



RESPIRATOR FITTING INSTRUCTIONS:

- 

STEP 1:
Hold the respirator in hand with the nose piece at your fingertips, allowing the headbands to hang freely below your hand.
- 

STEP 2:
Press the respirator firmly against your face with the nose piece on the bridge of your nose.
- 

STEP 3:
Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.
- 

STEP 4:
Using both hands, mold nose piece to the shape of your nose.
- 

STEP 5:
To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit.

Please carefully follow these fitting instructions during each use to achieve proper fit.

Lot Number/MFG Date:
Número de Lote/Fecha de Fabricación:

Expiration Date:
Fecha de Caducidad: