

SIZE:140*130*175mm

140 mm

130 mm

130 mm

175 mm

Makrite

9500V-N95



20
Particulate
Respirators



MAKRITE INDUSTRIES INC.
11F-5, No. 79 Sec. 1, Hsin Tai Wu Rd.
Hsichih, Taipei Hsien, Taiwan
Tel: (02)26962419



This respirator is approved only in the following configuration:

TC	Protection	Respirator	Cautions and Limitations ²
84A-XXXX	N95	9500V-N95	x ABCJMNOP

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply

2. Cautions and Limitations

- A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do not exceed maximum use concentrations established by regulatory standards.
J - Failure to properly use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations.
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.

SERVICE TIME LIMITATION

- Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

Makrite Industries, Inc. Phone No. in USA: 1-800-379-9929 www.makrite.com

Makrite

9500V-N95



heat welded straps



20
Particulate
Respirators

LIMITED WARRANTY
IMPORTANT NOTICE TO
PURCHASER:


This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Makrite will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Makrite products proved to be defectively manufactured. Except as provided above, Makrite shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

WARNINGS:


- Wearer must read fitting instructions and use limitations before use.
- This respirator does not protect against the risk of contracting disease or infection.
- Individuals highly sensitive to natural rubber latex may have an allergic reaction.
- Before use, must be fit tested either quantitatively or qualitatively by a qualified fit tester before using the respirator in a contaminated environment, as per the OSHA guidelines, 29 CFR 1910.134.

FITTING INSTRUCTIONS:

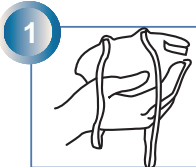

Step 1

Hold respirator in your hand as shown, with the nosepiece at your fingertips and the headbands hanging freely below your hand.

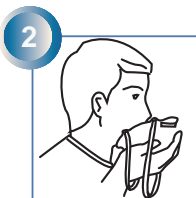


Step 2

Place respirator under your chin with the nosepiece up. Pull the top headband to top back of head.



Step 3

Pull the bottom headband overhead and place around neck below ears.

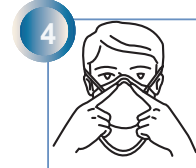

Step 4

Adjust headbands for comfortable fit. Push the nosepiece inward with your fingertips to conform to the shape of your nose.



Step 5

Cup both hands over the mask, checking for tightness. Inhale sharply, negative pressure should be felt in the mask. Re-adjust as necessary.

The respirator should be fit checked prior to each wearing. Place both hands completely over the respirator and exhale. If air leaks around nose or the edges, adjust the nosepiece and/or headbands until a good fit is achieved.

LOT NO.