

20
Disposable
Respirators

N95
Disposable
Respirators

**DIRECT
SUPPLY**
Taking Cost Out of Healthcare®
store.directsupply.com
1-800-634-7328
FQ942

**DIRECT
SUPPLY**
Taking Cost Out of Healthcare®
store.directsupply.com

1-800-634-7328
FQ942

DISPOSABLE N95 RESPIRATOR



20
Disposable
Respirators

SMALL

N95

NIOSH
Approved
Product

NIOSH National Institute for Occupational Safety and Health

DIRECT SUPPLY MANUFACTURING, INC.,
6635 West Champions Way, Milwaukee,
WI 53223, USA
Tel: +1414-491-6980

This respirator is approved only in the following configuration:

TC	Protection ¹	Respirator	Cautions and Limitations ²
84A-5463	N95	FQ942	2. ABCJMNOPS

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions apply.

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do not exceed maximum use concentrations established by regulatory standards.
J - Failure to properly use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used and maintained in accordance with MSHA, OSHA and other applicable regulations.
N - Never substitute, modify, add or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to user's instructions and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.
S - Special or critical User's Instructions and/or specific use limitation apply. Refer to User's Instructions before donning.

S - Special or Critical User's Instructions: The model FQ942 filtering facepiece respirator has been manufactured by Approval Holder Makrite for Company DIRECT SUPPLY MANUFACTURING, INC. under TC-84A-5463

SERVICE TIME LIMITATION

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
2. Leave contaminated area and remove respirator if breathing becomes difficult, proper fit cannot be obtained or the respirator is damaged, soiled or distorted. Do not use the same respirator again once it has been used in contaminated area.
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

LIMITED WARRANTY


IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, expressed or implied. There are no other warranties which extend beyond the description noted here. The physical standards and specifications of Direct Supply Manufacturing, Inc. will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Direct Supply Manufacturing, Inc. products proved to be defectively manufactured. Except as provided above, Direct Supply Manufacturing, Inc. shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

WARNINGS:

1. Wearer must read fitting instructions and limitations before use.
2. This respirator is not guaranteed to protect against the risk of contracting disease or infection.
3. Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training and fit testing.

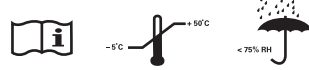

Made in China

**DIRECT
SUPPLY**
Taking Cost Out of Healthcare®
store.directsupply.com

1-800-634-7328
FQ942

DISPOSABLE N95 RESPIRATOR



20
Disposable
Respirators

SMALL

N95

NIOSH
Approved
Product

RESPIRATOR FITTING INSTRUCTION

- STEP 1**
Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.
- STEP 2**
Press the respirator firmly against your face with the nosepiece on the bridge of your nose.
- STEP 3**
Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.
- STEP 4**
Using both hands, mold nosepiece to the shape of your nose.
- STEP 5**
To test fit, cup both hands over the respirator while careful not to disturb placement, then inhale deeply. If air leaks around the edges, reposition the straps or adjust trap tension for better fit. Follow these fitting instructions during each use to achieve proper fit.

Caution:

The respirator should be fit checked prior to each use.

Product Description: Disposable N95 Respirator, NIOSH Approved, Small

Item Number: FQ942

Style Number: PM-DC20035

UOM: 20 ea/box

DIRECT SUPPLY MANUFACTURING, INC.,
Address: 6635 West Champions Way, Milwaukee, WI 53223, USA
Tel: +1414-491-6980
Manufactured by Dongguan Sure-On Plastic & Fiber Products Co., Ltd.
ZhuFo Ling Industrial Area, TangXia Town, DongGuan City, GuangDong, China
Tel: 86 769 8772 3170

Made in China

Executive Standard: 美国NIOSH 42 CFR Part 84 - 1995



LOT NO. / MFG Date:

Expiration Date: