

140 MM

130 MM

130 MM

140 MM

HF0804C



20 Pieces N95 Disposable Respirator

HF0804C



20 PIECES N95 Disposable Respirator

HF0804C




20 PIECES N95 Disposable Respirator



RESPIRATOR FITTING INSTRUCTIONS



STEP 1
· Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand

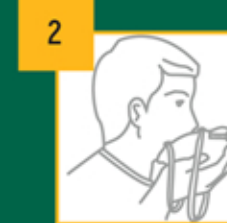

STEP 2
· Press the respirator firmly against your face with the nosepiece on the bridge of your nose



STEP 3
· Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.




STEP 4
· Using both hands, mold nosepiece to the shape of your nose



STEP 5
· To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

200 MM

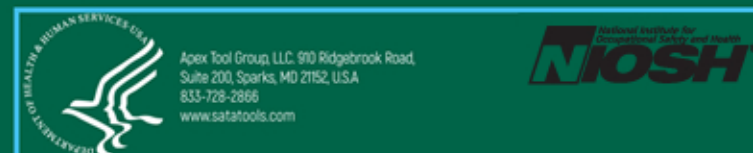

This respirator is approved only in the following configuration:

TC	Protection	Respirator	Cautions and Limitations ²
84A-5525	N95	HF0804C X	ABCI/MNOPS

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen

B - Not for use in atmospheres immediately dangerous to life or health

C - Do not exceed maximum use concentrations established by regulatory standards

D - Failure to properly use and maintain this product could result in injury or death

M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations

N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer

O - Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators

P - NIOSH does not evaluate respirators for use as surgical masks

S - Special or critical User's Instructions and/or specific use limitations apply. Refer to User's Instructions before donning

S - Special or Critical User's Instructions: The model HF0804C filtering facepiece respirator has been manufactured by Approval Holder Makrite for Apex Tool Group, LLC under TC-84A-5525

S - Special or critical User's Instructions: This product also offers nuisance level relief from organic vapors that are below the PEL (permissible exposure limit). Nuisance level refers to concentrations not exceeding OSHA PEL or other applicable government occupational exposure limits, whichever is lower.

SERVICE TIME LIMITATION

1. Do not use the same respirator for more than 8 hours, continuous or intermittent, in dry workplaces that could result in high filter loading.
2. Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once removed.
3. The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

LIMITED WARRANTY

IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, express or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Makrite will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Makrite products proved to be defectively manufactured. Except as provided above, Makrite shall not be liable or responsible for any loss, damages or liability, direct, indirect, incidental, special or consequential, arising out of sale, use, misuse or the inability to use products by the user.



Executive Standard: USA NIOSH 42 CFR Part 84 - 1995

LOT NO. / MFG Date:

Expiration Date: