## FLUIDSHIELD\* 3 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with Safety Seal and SO SOFT\* Lining

Orange

Small

- (FR) Masque chirurgical et respirateur à filtre à particules N95 FLUIDSHIELD\* 3 avec Safety Seal et doublure SO SOFT\*, orange, petit  
(DE) FLUIDSHIELD\* 3 N95 Partikelfilter- und OP-Maske mit Safety Seal und SO SOFT\*-Futter, orange, klein  
(BG) Респиратор с противопрахов филтър FLUIDSHIELD\* 3 №95 и хирургическа маска с подплата SO SOFT\*, оранжева, малък  
(ES) Mascarilla respiratoria y quirúrgica FLUIDSHIELD\* 3 con filtro de partículas N95, sello Safety Seal y forro SO SOFT\*, anaranjada, pequeña  
(CS) Respirátor FLUIDSHIELD\* 3 N95 s částicovým filtrem a chirurgickou maskou, těsněním Safety Seal a vložkou SO SOFT\*, oranžová, malá  
(DA) FLUIDSHIELD\* 3 kirurgmaske med N95-partikelfilter, Safety Seal og SO SOFT\* foring, orange, lille  
(ET) FLUIDSHIELD\* 3 N95 filter-respiraator ja kirurgiline mask Safety Seal ja SO SOFT\* voodriga, oranž, väike  
(EL) FLUIDSHIELD\* 3 N95 Αναπνευστική και χειρουργική μάσκα με φίλτρο σωματιδίων, μεμβράνη Safety Seal και επένδυση SO SOFT\*, πορτοκαλί, μικρό  
(IT) FLUIDSHIELD\* 3 N95 Respiratore antiparticolato e mascherina chirurgica con Safety Seal e rivestimento SO SOFT\*, arancione, piccolo  
(LV) Respirators ar FLUIDSHIELD\* 3 N95 sūkdaļu filtru un kirurgiskā maska, Safety Seal un SO SOFT\* odiri, oranža, mazs  
(LT) Chirurginė kaukė su FLUIDSHIELD\* 3 N95 dalelių filtru, Safety Seal ir SO SOFT\* pamušalu, oranžinė, maža  
(HU) FLUIDSHIELD\* 3 N95 részecszezsűrű respirátormasz és sebész maszk Safety Seal és SO SOFT\* béléssel, narancssárga, kis méret  
(NL) FLUIDSHIELD\* 3 N95 stoffiltermasker en operatiemasker met Safety Seal en SO SOFT\* beschermlaag, oranje, small

**210**  
(6 x 35)

MD

NON  
STERILE

30°C  
86°F  
-20°C  
-4°F

[halyardhealth.com/information](http://halyardhealth.com/information)

CE

\* The COLOR ORANGE is a  
Registered Trademark or  
Trademark of O&M Halyard, Inc.,  
or its affiliates.

REF

**46867 81**

O&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • [halyardhealth.com](http://halyardhealth.com)  
[EC] REP ArcRoyal Unlimited Company, Virginia Road, Kells; Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium  
Sponsored in Australia by O&M Halyard Australia Pty Ltd., 52 Alfred Street S, Milsons Point, NSW 2061  
製造販売元 O&M Halyard Japan合同会社 東京都港区芝公園2-6-3  
\*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-01 / 70215812

LOT

Product of USA

00000000000

UDI

Jul 31, 2020  
Thermal Case Label  
6" x 6.75"

## FLUIDSHIELD\* 3 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with Safety Seal and SO SOFT\* Lining

Orange

Small

- (NO) FLUIDSHIELD\* 3 N95 partikkelfiltrerende pustemaske og kirurgiskaske med Safety Seal og SO SOFT\* foring, oransje, liten
- (PL) Respirator z filtrem cząstek stałych FLUIDSHIELD\* 3 N95 i maska chirurgiczna oraz taśmą uszczelniającą Safety Seal i powłoką SO SOFT\*, pomarańczowa, rozmiar normalny
- (PT) Respirador com filtro de partículas FLUIDSHIELD\* 3 N95 e máscara cirúrgica com Safety Seal e com revestimento SO SOFT\*, cor-de-laranja, tamanho pequeno
- (RO) Mască chirurgicală și mască respiratorie cu filtru pentru particule FLUIDSHIELD\* 3 N95 cu Safety Seal și căptușeală SO SOFT\*, portocalii, mică
- (RU) Маска медицинская нестерильная одноразовая с пылепоглощающим фильтром FLUIDSHIELD\* 3 N95 и покрытием SO SOFT\*, оранжевая, малого размера
- (SK) Respirátor FLUIDSHIELD\* 3 N95 s časticovým filtrom a chirurgickou maskou, tesnením Safety Seal a vložkou SO SOFT\*, oranžový, malá
- (SL) FLUIDSHIELD\* 3 N95 respirator s filtrom za delce in kirurška maska in Safety Seal in SO SOFT\* oblogo, oranžen, majhen
- (FI) Hiukkassuodattimella varustettu FLUIDSHIELD\* 3 N95 -hengityssuojain ja leikkausmaski - Safety Seal ja SO SOFT\* -vuori, oranssi, pieni
- (SV) FLUIDSHIELD\* 3 N95 partikelfiltrerande andningsskydd och kirurgiskt munskydd med Safety Seal och SO SOFT\* foder, orange, litet
- (TR) FLUIDSHIELD\* 3 N95 Safety Seal SO SOFT\* Astarlı Toz Filtreli Solunum Cihazı ve Ameliyat Maskesi, turuncu, küçük
- (JA) FLUIDSHIELD\* 3 N95微粒子フィルターレスピレーター&サージカルマスク Safety Seal&SO SOFT\*ライナー、オレンジ、スモールサイズ
- (AR) كمامة FLUIDSHIELD\* 3 N95 مرشحة للجسيمات وقناع جراحي ومزودة Safety Seal وبطانة SO SOFT\*، برتقالية، صغيرة،
- (ZH) FLUIDSHIELD\* 3 N95微粒过滤呼吸器和手术口罩（含Safety Seal和SO SOFT\*内衬），橙色，小号

[halyardhealth.com/information](http://halyardhealth.com/information)

REF

# 4686781

O&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • [halyardhealth.com](http://halyardhealth.com)  
[EC] REP ArcRoyal Unlimited Company, Virginia Road, Kells; Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium  
Sponsored in Australia by O&M Halyard Australia Pty Ltd.; 52 Alfred Street S, Milsons Point, NSW 2061  
製造販売元 O&M Halyard 合同会社 東京都港区芝公園2-6-3  
\*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-01 / 70215812

LOT

Product of USA

000000000000

UDI

Jul 31, 2020  
Thermal Case Label  
6" x 6.75"


**FLUIDSHIELD\* 3**  
N95 PARTICULATE  
FILTER RESPIRATOR  
AND SURGICAL MASK

with SO SOFT\* Lining

Directions for Application

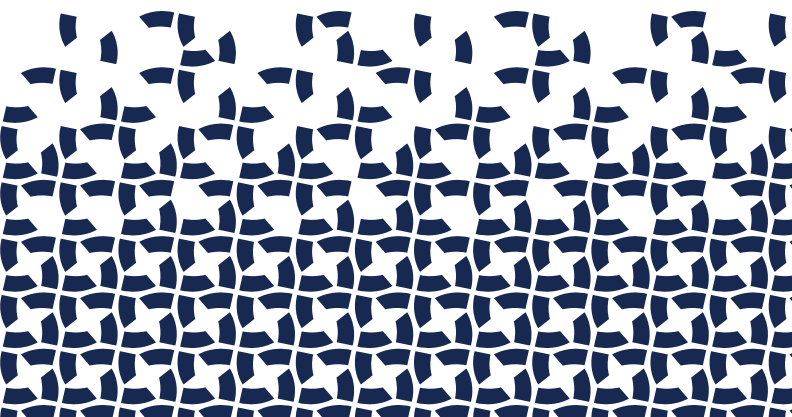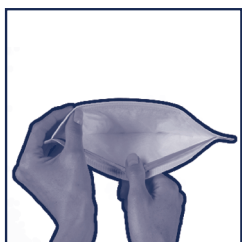
Fig. 1

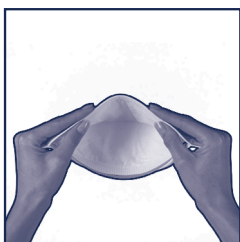
Fig. 2

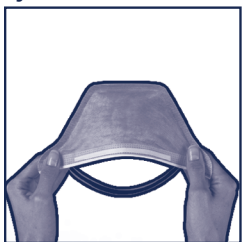
Fig. 3

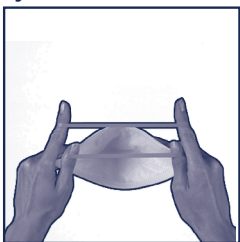
Fig. 4

Fig. 5

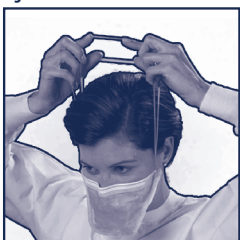
Fig. 6

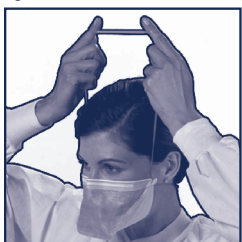
Fig. 7

Fig. 8

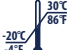
Fig. 9

Fig. 10

			
	 <a href="http://halyardhealth.com/information">halyardhealth.com/information</a>		

# HALYARD\* FLUIDSHIELD\* 3 N95 Particulate Filter Respirator and Surgical Mask with SO SOFT\* Lining

## Directions for Application

1. Separate the edges of the respirator to fully open it.
2. Slightly bend the nose wire to form a gentle curve.
3. Hold the respirator upside down to expose the two headbands.
4. Using your index fingers and thumbs, separate the two headbands.
5. While holding the headbands with your index fingers and thumbs, cup the respirator under your chin.
6. Pull the headbands up and over your head.
7. Release the lower headband from your thumbs and position it at the base of your neck.
8. Position the remaining headband on the crown of your head.
9. Conform the nose piece across the bridge of your nose by firmly pressing down with your fingers.
10. Continue to adjust the respirator and secure the edges until you feel you have achieved a good facial fit. Now, perform a Fit Check.

## Fitting Instructions:

- **Wear Orange side out.**
- Follow Directions for Application each time respirator is worn.
- Check facial seal to assure proper fit before entering any contaminated area.
  - Adjust nosepiece if there is air leakage around the nose.
  - Reposition headbands to secure facial seal around edges of respirator.
- Do not enter a contaminated area if a proper fit cannot be obtained. Follow all instructions and warnings on the use of this respirator, obtain a proper fit, and wear during all times of exposure. Failure to do so will reduce respirator effectiveness and user protection, and may result in sickness or death.

The HALYARD\* FLUIDSHIELD\* 3 N95 Particulate Filter Respirator and Surgical Masks (84A-7521, 84A-7523, 84A-7518, 84A-7520) are NIOSH approved as N95 particulate filter respirators. They provide 95% filtration efficiency of 0.3 micron particles. They meet the CDC Guidelines for TB exposure control.

## Intended Use:

The FLUIDSHIELD\* 3 N95 Particulate Filter Respirator and Surgical Mask is intended for use by operating room personnel and other health care workers to protect both patients and health care workers from transfer of microorganisms, blood and body fluids and airborne particulate materials.

## Use Instructions:

1. Before use, the wearer must be trained by the employer in proper respirator use in accordance with applicable safety and health standards.
2. If the respirator becomes damaged or breathing becomes difficult, leave the area. Discard and replace the respirator.

**WARNING: When properly worn, FLUIDSHIELD\* 3 N95 Particulate Filter Respirator and Surgical Mask reduces potential exposure to blood and body fluids as well as airborne microorganisms. This product does not completely eliminate the risk of contracting disease or infection.**

**REMINDER: Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing. Consult the facility's Respiratory Protection Program for further instruction.**

**These respirators are only approved in the following configurations:**

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-7521	N95	46727	ABCJMNP
84A-7523	N95	46767	
84A-7518	N95	46827	ABCJMNP
84A-7520	N95	46867	

## 1. PROTECTION:

**N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.**

## 2. CAUTIONS AND LIMITATIONS:

- A. Not for use in atmospheres containing less than 19.5% oxygen.
- B. Not for use in atmospheres immediately dangerous to life or health.
- C. Do Not exceed maximum use concentrations established by regulatory standards.
- J. Failure to properly use and maintain this product could result in injury or death.
- M. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O. Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators for use as surgical masks.

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