FLUIDSHIELD* 3 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with SO SOFT* Lining

Orange
Regular
Size

- (FR) Masque chirurgical et respirateur à filtre à particules N95 FLUIDSHIELD* 3 avec doublure SO SOFT*, orange, taille normale
(DE) FLUIDSHIELD* 3 N95 Partikelfilter- und OP-Maske mit SO SOFT*-Futter, orange, Standardgröße
(BG) Респиратор с противопрахов филтър FLUIDSHIELD* 3 №95 и хирургическа маска с подплата SO SOFT*, оранжева, стандартен
(ES) Mascarilla respiratoria y quirúrgica FLUIDSHIELD* 3 con filtro de partículas N95 y forro SO SOFT*, anaranjada, tamaño regular
(CS) Respirátor FLUIDSHIELD* 3 N95 s částicovým filtrem a chirurgickou maskou s vložkou SO SOFT*, oranžová, běžná velikost
(DA) FLUIDSHIELD* 3 kirurgmaske med N95-partikelfilter med SO SOFT* foring, orange, alm. størrelse
(ET) FLUIDSHIELD* 3 N95 filter-respiraator ja kirurgiline mask SO SOFT* voodriga, oranž, tavaline suurus
(EL) FLUIDSHIELD* 3 N95 Αναπνευστική και χειρουργική μάσκα με φίλτρο σωματιδίων και επένδυση SO SOFT*, πορτοκαλί, κανονικό μέγεθος
(IT) FLUIDSHIELD* 3 N95 Respiratore antiparticolato e mascherina chirurgica con rivestimento SO SOFT*, arancione, misura normale
(LV) Respirators ar FLUIDSHIELD* 3 N95 sīkdaļu filtru un ķirurģiskā maska un SO SOFT* oderī, oranžs, standarta izmērs
(LT) Chirurginė kaukė su FLUIDSHIELD* 3 N95 dalelių filtru ir SO SOFT* pamušalu, oranžinė, standartinis dydis
(HR) FLUIDSHIELD* 3 N95 respirator za filtriranje čestica i kirurška maska s podstavom SO SOFT*, narančasta, standardna veličina
(HU) FLUIDSHIELD* 3 N95 részecskeszűrő respirátormasz és sebész maszk SO SOFT* béleléssel, narancssárga, normál méret
(MT) Respiratur bil-filtru tal-partikuli N95 u maskra kirurgjika b'sigill ta' sigurtà FLUIDSHIELD* 3 b'inforra SO SOFT*, orangjo, daqs regolari
(NL) FLUIDSHIELD* 3 N95 stoffiltermasker en operatiemasker met SO SOFT* beschermelaag, oranje, standaardmaat

210
(6 x 35)

halyardhealth.com/information

*The COLOR ORANGE is a
Registered Trademark or
Trademark of O&M Halyard, Inc.,
or its affiliates.

REF

4672781

O&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • halyardhealth.com
[EC REP] ArcRoyal Unlimited Company, Virginia Road, Kells; Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium
Sponsored in Australia by O&M Halyard Australia Pty Ltd., 52 Alfred Street S, Milsons Point, NSW 2061
製造販売元 O&M Halyard Japan 合同会社 東京都港区芝公園2-6-3
*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-01 / 70215812

LOT

Product of USA

000000000000

UDI

May28, 2020
Thermal Case Label
6" x 6.75"

FLUIDSHIELD* 3 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with SO SOFT* Lining

Orange
Regular
Size

- (NO) FLUIDSHIELD* 3 N95 partikkelfiltrerende pustemaske og kirurgimaske med SO SOFT* foring, oransje, middels størrelse
(PL) Respirator z filtrem cząstek stałych FLUIDSHIELD* 3 N95 i maska chirurgiczna i powłoką SO SOFT*, pomarańczowa, rozmiar normalny
(PT) Respirador com filtro de partículas FLUIDSHIELD* 3 N95 e máscara cirúrgica com revestimento SO SOFT*, cor-de-laranja, tamanho regular
(RO) Mască chirurgicală și mască respiratorie cu filtru pentru particule FLUIDSHIELD* 3 N95 cu căptușeală SO SOFT*, portocalii, mărime normală
(RU) Маска медицинская нестерильная одноразовая с пылепоглощающим фильтром FLUIDSHIELD* 3 N95 и покрытием SO SOFT*, оранжевая, стандартного размера
(SK) Respirátor FLUIDSHIELD* 3 N95 s časticovým filtrom a chirurgickou maskou a vložkou SO SOFT*, oranžová, Bežná veľkosť
(SL) FLUIDSHIELD* 3 N95 respirator s filtrom za delce in kirurška maska in SO SOFT* oblogo, oranžen, običajna velikost
(SR) FLUIDSHIELD* 3 N95 mikročestični filterski respirator i hirurška maska sa SO SOFT* tkanjem, narandžasta, redovna veličina
(FI) Hiukkassuodattimella varustettu FLUIDSHIELD* 3 N95 -hengityssuojain ja leikkausmaski - SO SOFT* -vuori, oranssi, normaalikoko
(SV) FLUIDSHIELD* 3 N95 partikelfiltrerande andningsskydd och kirurgiskt munskydd med SO SOFT* foder, orange, normalstorlek
(TR) FLUIDSHIELD* 3 N95 SO SOFT* Astarlı Toz Filtreli Solunum Cihazı ve Ameliyat Maskesi, turuncu, normal boy
(JA) FLUIDSHIELD* 3 N95 微粒子フィルターレスピレーター&サージカルマスク
(KO) SO SOFT* 라이너, 오렌지, 레귤러 사이즈
(AR) FLUIDSHIELD* 3 N95 مرشحة للجسيمات وقناع جراحي ببطانة SO SOFT* حجم اعتيادي، برتقالية
(ZH) FLUIDSHIELD* N95 3 微粒过滤呼吸器和手术口罩 (含 SO SOFT* 内衬), 橙色, 标准大小

210
(6 x 35)

halyardhealth.com/information

REF

46727 81

O&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • halyardhealth.com
[EC REP] ArcRoyal Unlimited Company, Virginia Road, Kells; Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium
Sponsored in Australia by O&M Halyard Australia Pty Ltd., 52 Alfred Street S, Milsons Point, NSW 2061
製造販売元 O&M Halyard 合同会社 東京都港区芝公園2-6-3
*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-01 / 70215812

LOT

Product of USA

000000000000

UDI

May28, 2020
Thermal Case Label
6" x 6.75"


FLUIDSHIELD* 3
N95 PARTICULATE
FILTER RESPIRATOR
AND SURGICAL MASK

with SO SOFT* Lining

Directions for Application

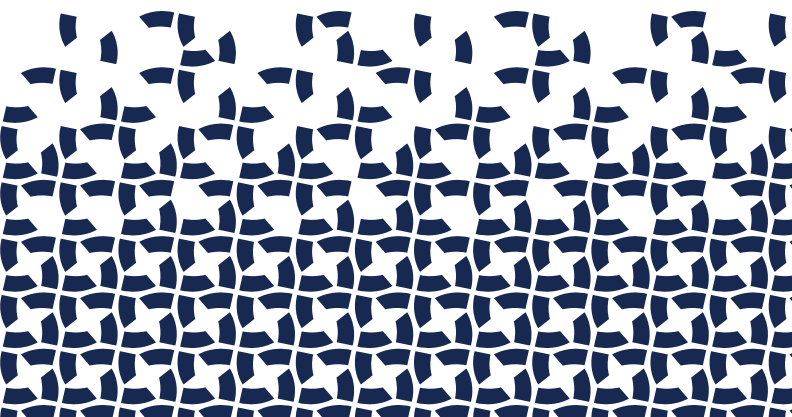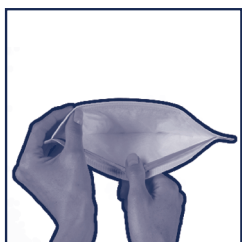
Fig. 1

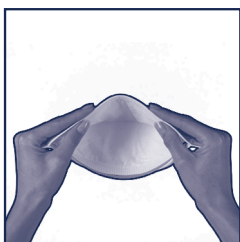
Fig. 2

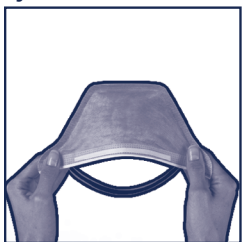
Fig. 3

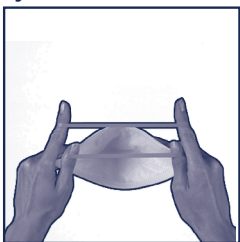
Fig. 4

Fig. 5

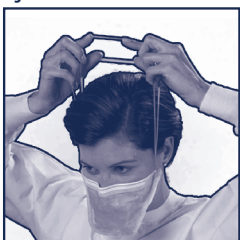
Fig. 6

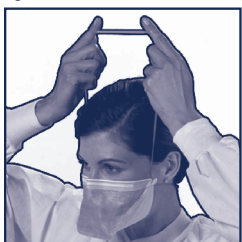
Fig. 7

Fig. 8

Fig. 9

Fig. 10

	halyardhealth.com/information		

HALYARD* FLUIDSHIELD* 3 N95 Particulate Filter Respirator and Surgical Mask with SO SOFT* Lining

Directions for Application

1. Separate the edges of the respirator to fully open it.
2. Slightly bend the nose wire to form a gentle curve.
3. Hold the respirator upside down to expose the two headbands.
4. Using your index fingers and thumbs, separate the two headbands.
5. While holding the headbands with your index fingers and thumbs, cup the respirator under your chin.
6. Pull the headbands up and over your head.
7. Release the lower headband from your thumbs and position it at the base of your neck.
8. Position the remaining headband on the crown of your head.
9. Conform the nose piece across the bridge of your nose by firmly pressing down with your fingers.
10. Continue to adjust the respirator and secure the edges until you feel you have achieved a good facial fit. Now, perform a Fit Check.

Fitting Instructions:

- **Wear Orange side out.**
- Follow Directions for Application each time respirator is worn.
- Check facial seal to assure proper fit before entering any contaminated area.
 - Adjust nosepiece if there is air leakage around the nose.
 - Reposition headbands to secure facial seal around edges of respirator.
- Do not enter a contaminated area if a proper fit cannot be obtained. Follow all instructions and warnings on the use of this respirator, obtain a proper fit, and wear during all times of exposure. Failure to do so will reduce respirator effectiveness and user protection, and may result in sickness or death.

The HALYARD* FLUIDSHIELD* 3 N95 Particulate Filter Respirator and Surgical Masks (84A-7521, 84A-7523, 84A-7518, 84A-7520) are NIOSH approved as N95 particulate filter respirators. They provide 95% filtration efficiency of 0.3 micron particles. They meet the CDC Guidelines for TB exposure control.

Intended Use:

The FLUIDSHIELD* 3 N95 Particulate Filter Respirator and Surgical Mask is intended for use by operating room personnel and other health care workers to protect both patients and health care workers from transfer of microorganisms, blood and body fluids and airborne particulate materials.

Use Instructions:

1. Before use, the wearer must be trained by the employer in proper respirator use in accordance with applicable safety and health standards.
2. If the respirator becomes damaged or breathing becomes difficult, leave the area. Discard and replace the respirator.

WARNING: When properly worn, FLUIDSHIELD* 3 N95 Particulate Filter Respirator and Surgical Mask reduces potential exposure to blood and body fluids as well as airborne microorganisms. This product does not completely eliminate the risk of contracting disease or infection.

REMINDER: Before occupational use of this respirator, a written respiratory protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing. Consult the facility's Respiratory Protection Program for further instruction.

These respirators are only approved in the following configurations:

TC-	Protection ¹	Respirator	Cautions and Limitations ²
84A-7521	N95	46727	ABCJMNP
84A-7523	N95	46767	
84A-7518	N95	46827	ABCJMNP
84A-7520	N95	46867	

1. PROTECTION:

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.

2. CAUTIONS AND LIMITATIONS:

- A. Not for use in atmospheres containing less than 19.5% oxygen.
- B. Not for use in atmospheres immediately dangerous to life or health.
- C. Do Not exceed maximum use concentrations established by regulatory standards.
- J. Failure to properly use and maintain this product could result in injury or death.
- M. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O. Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators for use as surgical masks.

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