# FLUIDSHIELD\* 2 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with SO SOFT\* Lining

White

Small

- (FR) Masque chirurgical et respirateur à filtre à particules N95 FLUIDSHIELD\* 2 avec doublure SO SOFT\*, blanc, petit  
(DE) FLUIDSHIELD\* 2 N95 Partikelfilter- und OP-Maske mit SO SOFT\*-Futter, weiss, klein  
(BG) Респиратор с противопрахов филтър FLUIDSHIELD\* 2 N95 и хирургическа маска с подплата SO SOFT\*, бяла, малък  
(ES) Mascarilla respiratoria y quirúrgica con filtro de partículas FLUIDSHIELD\* 2 N95 y forro SO SOFT\*, blanca, pequeña  
(CS) Respirátor FLUIDSHIELD\* 2 N95 s částicovým filtrem a chirurgickou maskou s vložkou SO SOFT\*, bílý, malý  
(DA) Kirurgimaske med FLUIDSHIELD\* 2 N95-partikelfilter med SO SOFT\* foring, hvid, lille  
(ET) FLUIDSHIELD\* 2 N95 filter-respiraator ja kirurgiline mask SO SOFT\* voodriga, valge, väike  
(EL) FLUIDSHIELD\* 2 N95 Αναπνευστική και χειρουργική μάσκα με φίλτρο σωματιδίων και επένδυση SO SOFT\*, λευκή, μικρό μέγεθος  
(IT) FLUIDSHIELD\* 2 N95 Respiratore antiparticolato e mascherina chirurgica con rivestimento SO SOFT\*, bianco, piccolo  
(LV) Respirators ar FLUIDSHIELD\* 2 N95 sikdāju filtru un ķirurģiskā maska ar SO SOFT\* oderi, balts, mazs  
(LT) Chirurginė kaukė su FLUIDSHIELD\* 2 N95 dalelių filtru ir SO SOFT\* pamušalu, balta, maža  
(HU) FLUIDSHIELD\* 2 N95 részecskeszűrő respirátormaszk és sebészeti maszk SO SOFT\* béléssel, fehér, kis méret  
(NL) FLUIDSHIELD\* 2 N95 stoffiltermasker en operatiemasker met SO SOFT\* bescherm laag, wit, small

**300**  
(6 x 50)

[halyardhealth.com/information](http://halyardhealth.com/information)

REF

**62355 81**

Q&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • [halyardhealth.com](http://halyardhealth.com)  
[EC] REP ArcRoyal Unlimited Company, Virginia Road, Kells, Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium  
Sponsored in Australia by O&M Halyard Australia Pty Ltd., Tenancy 1, Level 1, Trinit 3, 39 Delhi Road, North Ryde NSW 2113  
製造販売元 O&M Halyard Japan 合同会社 東京都港区芝公園2-6-3 芝公園フロントタワー12階  
\*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-02 / 70216790

LOT

Product of USA

UDI

000000000000

Jul 31, 2020  
Thermal Case Label  
6" x 6.75"

## FLUIDSHIELD\* 2 N95 PARTICULATE FILTER RESPIRATOR AND SURGICAL MASK

with SO SOFT\* Lining

White

Small

- (NO) FLUIDSHIELD\* 2 N95 Partikkelfiltrerende pustemaske og kirurgimaske med SO SOFT\* foring, hvit, liten  
(PL) Respirator z filtrem cząstek stałych FLUIDSHIELD\* 2 N95 i maska chirurgiczna z powłoką SO SOFT\*, biała, mała  
(PT) Respirador com filtro de partículas FLUIDSHIELD\* 2 N95 e máscara cirúrgica com revestimento SO SOFT\*, branco, tamanho pequeno  
(RO) Mască chirurgicală și mască respiratorie cu filtru pentru particule FLUIDSHIELD\* 2 N95 cu căptușeală SO SOFT\*, alb, mic  
(RU) Маска медицинская нестерильная одноразовая с пылепоглощающим фильтром FLUIDSHIELD\* 2 N95 и покрытием SO SOFT\*, белая, малого размера  
(SK) Dýchací prístroj FLUIDSHIELD\* 2 N95 s časticovým filtrom a chirurgickou maskou s vložkou SO SOFT\*, biely, malý  
(SL) FLUIDSHIELD\* 2 N95 respirator s filtrom za delce in kirurška maska s SO SOFT\* oblogo, bel, majhen  
(FI) Hiukkassuodattimella varustettu FLUIDSHIELD\* 2 N95 -hengityssuojain ja leikkausmaski - SO SOFT\*-vuori, valkoinen, pieni  
(SV) FLUIDSHIELD\* 2 N95 partikelfiltrerande andningsskydd och kirurgiskt munskydd, med SO SOFT\* foder, vitt, litet  
(TR) FLUIDSHIELD\* 2 N95 SO SOFT\* Astarlı Toz Filtreli Solunum Cihazı ve Ameliyat Maskesi, beyaz, küçük  
(JA) FLUIDSHIELD\* 2 N95微粒子フィルターレスピレーター&サージカルマスク SO SOFT\*ライナー、白、スモールサイズ  
(AR) كمامة N95 مرشحة للجسيمات وقناع جراحي ببطانة So Soft\*، بيضاء، صغيرة  
(ZH) FLUIDSHIELD\* 2 N95微粒过滤呼吸器和手术口罩 (含SO SOFT\*内衬)，白色，小号

**300**  
(6 x 50)

[halyardhealth.com/information](http://halyardhealth.com/information)

REF

**62355 81**

O&M Halyard, Inc., 9120 Lockwood Blvd., Mechanicsville, VA 23116; In USA, please call 1-844-425-9273 • [halyardhealth.com](http://halyardhealth.com)  
[EC] [REP] ArcRoyal Unlimited Company, Virginia Road, Kells, Co Meath, Ireland • O&M Halyard Belgium BVBA, Berkenlaan 8B, 1831 Machelen (Brab.), Belgium  
Sponsored in Australia by O&M Halyard Australia Pty Ltd., Tenancy 1, Level 1, Trinitri 3, 39 Delhi Road, North Ryde NSW 2113  
製造販売元 O&M Halyard Japan合同会社 東京都港区芝公園2-6-3 芝公園フロントタワー12階  
\*Registered Trademark or Trademark of O&M Halyard, Inc., or its affiliates. ©2018. All rights reserved. 22-HM-001-0-02 / 70216790

LOT

Product of USA

000000000000

UDI

Jul 31, 2020  
Thermal Case Label  
6" x 6.75"


**FLUIDSHIELD\* 2**  
**N95 PARTICULATE**  
**FILTER RESPIRATOR**  
**AND SURGICAL MASK**

with SO SOFT\* Lining

**Directions for Application**

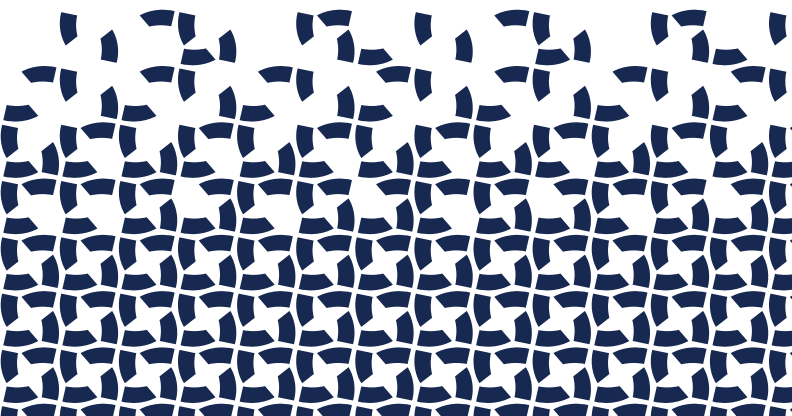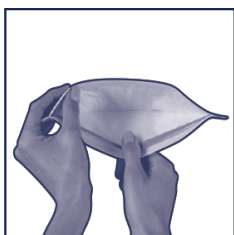
Fig. 1

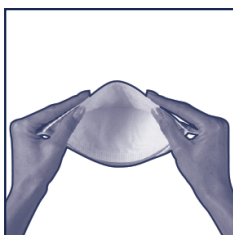
Fig. 2

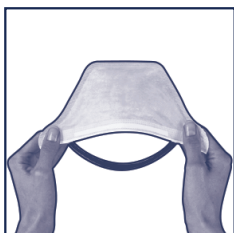
Fig. 3

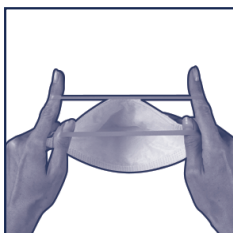
Fig. 4

Fig. 5

Fig. 6

Fig. 7

Fig. 8

Fig. 9

Fig. 10

	<a href="http://halyardhealth.com/information">halyardhealth.com/information</a>		

# **HALYARD\* FLUIDSHIELD\* 2** **N95 Particulate Filter Respirator and Surgical Mask** **with SO SOFT\* Lining**

## **Directions for Application**

1. Separate the edges of the respirator to fully open it.
2. Slightly bend the nose wire to form a gentle curve.
3. Hold the respirator upside down to expose the two headbands.
4. Using your index fingers and thumbs, separate the two headbands.
5. While holding the headbands with your index fingers and thumbs, cup the respirator under your chin.
6. Pull the headbands up and over your head.
7. Release the lower headband from your thumbs and position it at the base of your neck.
8. Position the remaining headband on the crown of your head.
9. Conform the nose piece across the bridge of your nose by firmly pressing down with your fingers.
10. Continue to adjust the respirator and secure the edges until you feel you have achieved a good facial fit. Now, perform a Fit Check.

## **Fitting Instructions:**

- Follow Directions for Application each time respirator is worn.
- Check facial seal to assure proper fit before entering any contaminated area.
  - Adjust nosepiece if there is air leakage around the nose.
  - Reposition headbands to secure facial seal around edges of respirator.
- Do not enter a contaminated area if a proper fit cannot be obtained. Follow all instructions and warnings on the use of this respirator, obtain a proper fit, and wear during all times of exposure. Failure to do so will reduce respirator effectiveness and user protection, and may result in sickness or death.

The HALYARD\* FLUIDSHIELD\* 2 N95 Particulate Filter Respirator and Surgical Masks (84A-7524, 84A-7525) are NIOSH approved as N95 particulate filter respirators. They provide 95% filtration efficiency of 0.3 micron particles. They meet the CDC Guidelines for TB exposure control.

## **Intended Use:**

The FLUIDSHIELD\* 2 N95 Particulate Filter Respirator and Surgical Mask is intended for use by operating room personnel and other health care workers to protect both patients and health care workers from transfer of microorganisms, blood and body fluids and airborne particulate materials.

## **Use Instructions:**

1. Before use, the wearer must be trained by the employer in proper respirator use in accordance with applicable safety and health standards.
2. If the respirator becomes damaged or breathing becomes difficult, leave the area. Discard and replace the respirator.

**WARNING: When properly worn, FLUIDSHIELD\* 2 N95 Particulate Filter Respirator and Surgical Mask reduces potential exposure to blood and body fluids as well as airborne microorganisms. This product does not completely eliminate the risk of contracting disease or infection.**

**REMINDER: Before occupational use of this respirator, a written protection program must be implemented meeting all the local government requirements. In the United States, employers must comply with OSHA 29 CFR 1910.134 which includes medical evaluation, training, and fit testing. Consult the facility's Respiratory Protection Program for further instruction.**

**These respirators are only approved in the following configurations:**

TC-	Protection <sup>1</sup>	Respirator	Cautions and Limitations <sup>2</sup>
84A-7524	N95	62355	ABCJMNOP
84A-7525	N95	62126	ABCJMNOP

## **1. PROTECTION:**

**N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil; time use restrictions may apply.**

## **2. CAUTIONS AND LIMITATIONS:**

- A. Not for use in atmospheres containing less than 19.5% oxygen.
- B. Not for use in atmospheres immediately dangerous to life or health.
- C. Do Not exceed maximum use concentrations established by regulatory standards.
- J. Failure to properly use and maintain this product could result in injury or death.
- M. All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA, and other applicable regulations.
- N. Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
- O. Refer to user's instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
- P. NIOSH does not evaluate respirators for use as surgical masks.

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