

尺寸：140X130X180MM
颜色：CMYK

40
N95 Disposable Respirators

GN8019-N95



GN8019-N95

NIOSH
Approved
Product

N95 Disposable Respirator

95% Efficient Protection Against
Non-Oil Particulates



40
Disposable
Respirators

LOT NO. / MFG Date:

Expiration Date:

Xuan Cheng Guang Neng Non-woven Co., Ltd.
58 Guangming Avenue, Economic
Development Zone, Xuanzhou District,
Xuancheng City, Anhui Province, China.
Tel: 0086 0563-3497599

TC ²	Protection ¹	Respirator	Cautions and Limitations ²
84A-7638	N95	GN8019-N95	ABCJMNOPS

1. Protection

N95-Particulate Filter (95% filter efficiency level) effective against particulate aerosols free of oil-time use restrictions may apply

2. Cautions and Limitations

A - Not for use in atmospheres containing less than 19.5% oxygen.
B - Not for use in atmospheres immediately dangerous to life or health.
C - Do Not exceed maximum use concentrations established by regulatory standards
J - Failure to properly use and maintain this product could result in injury or death.
M - All approved respirators shall be selected, fitted, used, and maintained in accordance with MSHA, OSHA and other applicable regulations
N - Never substitute, modify, add, or omit parts. Use only exact replacement parts in the configuration as specified by the manufacturer.
O - Refer to users instructions, and/or maintenance manuals for information on use and maintenance of these respirators.
P - NIOSH does not evaluate respirators for use as surgical masks.
S - Special or critical User's Instructions and/or specific use limitation apply.
Refer to User's Instructions before donning

S - Special or Critical User's Instructions: The model GN801-N95 filtering facepiece respirator has been manufactured by Approval Holder Makrite for Company XUAN CHENG GUANG NENG NON-WOVEN CO., LTD under TC-84A-7638.

SERVICE TIME LIMITATION


- Do not use the same respirator for more than 8 hours, continuous or intermittent, in dirty workplaces that could result in high filter loading.
- Leave contaminated area and remove respirator if breathing becomes difficult or proper fit cannot be obtained or the respirator is damaged, soiled, or distorted. Do not use the same respirator again once it has been used in contaminated area.
- The maximum service time of this respirator is five (5) consecutive calendar days (including days of non-use), beginning from the first day of use.

LIMITED WARRANTY IMPORTANT NOTICE TO PURCHASER:

This limited warranty is made in lieu of the warranties of merchantability, fitness for particular purposes and all other warranties, expressed or implied. There are no other warranties which extend beyond the description on the face hereof. The physical standards and specifications of Guangneng will be met by products sold.

Exclusive Remedies: damages for the breach of this limited warranty are limited to the replacement of such quantity of Guangneng, products proved to be defectively manufactured. Except as provided above, Guangneng shall not be liable or responsible for any loss, damage, or liability, direct, indirect, incidental, special, or consequential, arising out of sale, use, or misuse, or the inability to use products by the user.

Made in China Non-medical use

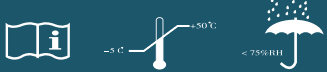

GN8019-N95

NIOSH
Approved
Product

N95 Disposable Respirator

95% Efficient Protection Against
Non-Oil Particulates



40
Disposable
Respirators

RESPIRATOR FITTING INSTRUCTION

- STEP 1**
Hold the respirator in hand with the nosepiece at your fingertips, allowing the headbands to hang freely below your hand.
- STEP 2**
Press the respirator firmly against your face with the nosepiece on the bridge of your nose.
- STEP 3**
Stretch and position the top band high on the back of the head. Stretch the bottom band over the head and position below your ears.
- STEP 4**
Using both hands, mold nosepiece to the shape of your nose.
- STEP 5**
To test fit: a) Cup both hands over the respirator being careful not to disturb position, and b) inhale vigorously. If air leaks around the edges, reposition the straps or adjust strap tension for better fit. Please carefully follow these fitting instructions during each use to achieve proper fit.

Xuan Cheng Guang Neng Non-woven Co., Ltd.
Address: 58 Guangming Avenue, Economic Development Zone, Xuanzhou District, Xuancheng City, Anhui Province, China.
Tel: 0086 0563-3497599
Manufactured by Dongguan Sure_On Plastic & Fiber Products Co., Ltd.
ZhuFo Ling Industrial Area, TangXia Town, DongGuan City, GuangDong, China
Tel: 86 769 8772 3170

Made in China